
WORLD INTELLECTUAL PROPERTY ORGANIZATION

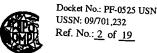

INTERNATIONAL APPLICATION PUBLISHED UNDER THE PATENT COOPERATION TREATY (PCT)

(51) International Patent Classification 6: C12Q 1/68, G06F 15/00 A1	(11) International Publication Number: WO 95/2068] (43) International Publication Date: 3 August 1995 (03.08.95
(21) International Application Number: (21) PCT/US95/0116	
(22) International Filing Date: 27 January 1995 (27:01:95	EE, FI, GE, HU, JP, KG, KP, KR, KZ, LK, LR, LT, LV MD, MG, MN, MX, NO, NZ, PL, RO, RU, SI, SK, TJ, TT UA, UZ, VN, European patent (AT, BE, CH, DE, DK, ES FR, GB, GR, IE, IT, LU, MC, NL, PT, SE), OAPI paten
(30) Priority Data: 08/187,530 27 January 1994 (27.01.94) U 08/282,955 29 July 1994 (29.07.94) U	TG). ARIPO natent (KF MW SD SZ)
(71) Applicant: INCYTE PHARMACEUTICALS, INC. [US/US] 3330 Hillview Avenue, Palo Alto, CA 94304 (US).	With international search report.
72) Inventors: SEILHAMER, Jeffrey, J.; 12555 La Cresta, Lo Altos Hills, CA 94022 (US). SCOTT, Randal, W.; 1314 Sun-Mor, Mountain View, CA 94040 (US).	of the thinking a commentation and a
(74) Agents: CAGE, Kenneth, L. et al.; Willian Brinks Hofer Gilson & Lione, 2000 K Street, N.W., Suite 200, Washington, DC 20006-1809 (US).	tion that they decided appears to
	 (1) 自然のできる。 View Segment (2) は、 (2) を表している。 View Segment (2) には、 (2) を必要できる。 (2) を表している。
	The first of the ending of the property of the contract of the
en de la companya de La companya de la co	e far the second of the second

(57) Abstract

A method and system for quantifying the relative abundance of gene transcripts in a biological specimen. One embodiment of the method generates high-throughput sequence-specific analysis of multiple RNAs or their corresponding cDNAs (gene transcript imaging analysis). Another embodiment of the method produces a gene transcript imaging analysis by the use of high-throughput cDNA sequence analysis. Another embodiment of the method produces a gene transcript imaging analysis by the use of high-unoughput come sequence analysis. In addition, the gene transcript imaging can be used to detect or diagnose a particular biological state, disease, or condition which is correlated to the relative abundance of gene transcripts in a given cell or population of cells. The invention provides a method for comparing the gene transcript image analysis from two or more different biological specimens in order to distinguish between the two specimens and identify one or more genes which are differentially expressed between the two specimens.

- The high standing on his and control (1900) is the set of all and only to act of the control of the control
- (a) The Control of the second blues of the community of the open page of the community of the community of the community of the control of
- (4) The Common of the Commo
- が、 Andrew Andr

FOR THE PURPOSES OF INFORMATION ONLY

Codes used to identify States party to the PCT on the front pages of pamphlets publishing international applications under the PCT.

ΑT	Austria	GB	United Kingdom	MIR	Mauritania
U	Australia	GE	Georgia	MW	Malawi
3B	Barbados	GN	Guinea	NE	Niger
BE	Belgium	GR	Greece	NL	Netherlands
BF	Burkina Faso	HU	Hungary	NO	Norway
BG	Bulgaria	IE	Ireland	NZ	New Zealand
BJ	Benin	IT	Italy	PL	Poland
BR	Brazil	JP	Japan	PT	Portugal
BY	Belarus	KE	Kenya	RO	Romania
CA	Canada	KG	Kyrgystan	RU	Russian Federation
CF	Central African Republic	KP	Democratic People's Republic	SD	Sudan
CG	Congo		of Korea	SE	Sweden
CH	Switzerland	KR	Republic of Korea	SI	Slovenia
CI	Côte d'Ivoire	KZ	Kazakhstan	SK	Slovakia
CM	Carneroon	LI	Liechtenstein	SN	Senegal
CN	China	LK	. Sri Lanka	TD	Chad
CS	Czechoslovakia	LU	Luxembourg	TG	Togo
CZ	Czech Republic .	LV	Larvia	TJ	Tajikistan
DE	Germany	MC	Monaco	TT	Trinidad and Tobago
DK	Denmark	MD	Republic of Moldova	UA	Ukraine
ES	Spain	MG	Madagascar	US	United States of America
FI	Finland	ML	Mali	UZ	Uzbekistan
FR	France	MN	Mongolia	VN	Viet Nam
GA ·	Gabon				

established the same of their same of the same of the same and the same of the

COMPARATIVE GENE TRANSCRIPT ANALYSIS

Proceed and the pale Field Of Invention and a continued

The present invention is in the field of molecular biology and computer science; more particularly, the present invention describes methods of analyzing gene transcripts and diagnosing the genetic expression of cells and tissue.

Description of the invention and the invention

Until very recently, the history of molecular biology
has been written one generatiantime. Scientists have
observed the cell's physical changes, isolated mixtures
from the cell or its milieu, purified proteins, sequenced
proteins and therefrom constructed probes to look for the
corresponding gene.

Recently, different nations have set up massive projects to sequence the billions of bases in the human genome. These projects typically begin with dividing the genome into large portions of chromosomes and then determining the sequences of these pieces, which are then analyzed for identity with known proteins or portions thereof, known as motifs. Unfortunately, the majority of genomic DNA does not encode proteins and though it is postulated to have some effect on the cell's ability to make protein, its relevance to medical applications is not understood at this time.

A third methodology involves sequencing only the transcripts encoding the cellular machinery actively involved in making protein, namely the mRNA. The advantage is that the cell has already edited out all the non-coding DNA, and it is relatively easy to identify the protein-coding portion of the RNA. The utility of this approach was not immediately obvious to genomic researchers. In fact, when cDNA sequencing was initially proposed, the method was roundly denounced by those committed to genomic sequencing. For example, the head of the U.S. Human Genome project discounted CDNA sequencing as not valuable and refused to approve funding of projects.

In this disclosure, we teach methods for analyzing DNA, including cDNA libraries. Based on our analyses and

research, we see each individual gene product as a "pixel" of information, which relates to the expression of that, and only that, gene. We teach herein, methods whereby the individual "pixels" of gene expression information can be combined into a single gene transcript "image," in which each of the individual genes can be visualized simultaneously and allowing relationships between the gene pixels to be easily visualized and understood.

subtraction. Electronic subtraction will enable the gene researcher to turn a single image into a moving picture, one which describes the temporality or dynamics of gene expression, at the level of a cell or a whole tissue. It is that sense of "motion" of cellular machinery on the scale of a cell or organ which constitutes the new invention herein. This constitutes a new view into the process of living cell physiology and one which holds great promise to unveil and discover new therapeutic and diagnostic approaches in medicine.

We teach another method which we call "electronic northern," which tracks the expression of a single gene across many types of cells and tissues.

20

Nucleic acids (DNA and RNA) carry within their sequence the hereditary information and are therefore the prime molecules of life. Nucleic acids are found in all living organisms including bacteria, fungi, viruses, plants and animals. It is of interest to determine the relative abundance of different discrete nucleic acids in different cells, tissues and organisms over time under various conditions, treatments and regimes.

All dividing cells in the human body contain the same set of 23 pairs of chromosomes. It is estimated that these autosomal and sex chromosomes encode approximately 100,000 genes. The differences among different types of cells are believed to reflect the differential expression of the 100,000 or so genes. Fundamental questions of biology could be answered by understanding which genes are transcribed and knowing the relative abundance of transcripts in different cells.

Previously, the art has only provided for the analysis of a few known genes at a time by standard molecular biology techniques such as PCR, northern blot analysis, or other types of DNA probe analysis such as in situ , ... 5 hybridization. Each of these methods allows one to analyze the transcription of only known genes and/or small numbers of genes at a time. Nucl. Acids Res. 19, 7097-7104 (1991); Nucl. Acids Res. 18, 4833-42 (1990); Nucl. Acids Res. 18, 4 2789-92 (1989); European J. Neuroscience 2, 1063-1073 . . 10 (1990); Analytical Biochem. 187,0 364-73 (1990); Genet. 187,0 364-73 Annals Techn. Appl. 7, 64-70 (1990); GATA 8(4), 129-33 (1991); Proc. Natl. Acad. Sci. USA 85, 1696-1700 (1988); Nucl. Acids Res. 19, 1954 (1991); Proc. Natl. Acad. Sci. USA 88, 1943-47 (1991); Nucl. Acids Res. 19, 6123-27 (1991) 15 (1991); Proc. Natl. Acad. Sci. USA 85, 5738-42 (1988); (1988) Studies of the number and types of genes whose transcription is induced or otherwise regulated during cell processes such as activation, differentiation, aging, viral 20 transformation, morphogenesis, and mitosis have been pursued for many years, using a variety of methodologies. One of the earliest methods was to isolate and analyze levels of the proteins in a cell, tissue, organ system, for even organisms both before and after the process of interest. One method of analyzing multiple proteins in a 25 sample is using 2-dimensional gel electrophoresis, wherein proteins can be, in principle, identified and quantified as individual bands, and ultimately reduced to a discrete At present, 2-dimensional analysis only resolves approximately 15% of the proteins. In order to positively analyze those bands which are resolved, each band must be excised from the membrane and subjected to protein sequence analysis using Edman degradation. Unfortunately, most of the bands were present in quantities too small to obtain a reliable sequence, and many of those bands contained more than one discrete protein. An additional difficulty is that many of the proteins were blocked at the amino-terminus, further complicating the sequencing

process.

Analyzing differentiation at the gene transcription as level has overcome many of these disadvantages and drawbacks, since the power of recombinant DNA technology allows amplification of signals containing very small 5 amounts of material. The most common method, called "hybridization subtraction," involves isolation of mRNA from the biological specimen before (B) and after (A) the developmental process of interest, transcribing one set of mRNA into cDNA, subtracting specimen B from specimen A 10 (mRNA from cDNA) by hybridization, and constructing a cDNA library from the non-hybridizing mRNA fraction. Many different groups have used this strategy successfully; and a variety of procedures have been published and improved upon using this same basic scheme. Nucl. Acids Res. 19, 15 7097-7104 (1991); Nuclea Acids Res. 18, 4833-420(1990); • Nucl. Acids Res. <u>18</u>, 2789-92 (1989); European J. Westurne Neuroscience 2, 1063-1073 (1990); Analytical Biochem. 187; 364-73 (1990); Genet. Annals Techn. Appl. 7, 64-70 (1990); GATA 8(4), 129-33 (1991); Proc. Natl. Acad. Sci. USA 85, 44 1696-1700 (1988); Nucl. Acids Res. 19, 1954 (1991); Proc. 20 Natl. Acad. Sci. USA 88, 1943-47 (1991); Nucl. Acids Res. 19, 6123-27 (1991); Proc. Natl. Acad. Sci. USA 85, 5738-42 (1988); Nucl. Acids Res. 16, 10937 (1988).

Although each of these techniques have particular strengths and weaknesses, there are still some limitations and undesirable aspects of these methods: First, the time and effort required to construct such libraries is quite large. Typically, a trained molecular biologist might expect construction and characterization of such a library 30 to require 3 to 6 months, depending on the level of skill, experience, and luck. Second, the resulting subtraction libraries are typically inferior to the libraries constructed by standard methodology. A typical conventional cDNA library should have a clone complexity of at least 106 clones, and an average insert size of 1-3 kB. In contrast, subtracted libraries can have complexities of 10^2 or 10^3 and average insert sizes of 0.2 kB. there can be a significant loss of clone and sequence information associated with such libraries. Third, this

approach allows the researcher to capture only the genes induced in specimen A relative to specimen B, not vice-versa, nor does it easily allow comparison to a third specimen of interest (C). Fourth, this approach requires very large amounts (hundreds of micrograms) of "driver" mRNA (specimen B), which significantly limits the number and type of subtractions that are possible since many tissues and cells are very difficult to obtain in large quantities.

upon the physical properties of DNA:DNA or RNA:DNA hybridization. The ability of a given sequence to find a hybridization match is dependent on its unique CoT value.

The CoT value is a function of the number of copies:

(concentration) of the particular sequence, multiplied by
the time of hybridization. It follows that for sequences
which are abundant, hybridization events will occur very
rapidly (low CoT value), while rare sequences will form the
duplexes at very high CoT values. CoT values which allow

such rare sequences to form duplexes and therefore be
effectively selected are difficult to achieve in a
convenient time frame. Therefore, hybridization
subtraction is simply not a useful technique with which to
study relative levels of rare mRNA species. Sixth, this
problem is further complicated by the fact that duplex

formation is also dependent on the nucleotide base composition for a given sequence. Those sequences rich in G + C form stronger duplexes than those with high contents of A + T. Therefore, the former sequences will tend to be

removed selectively by hybridization subtraction. Seventh, it is possible that hybridization between nonexact matches can occur. When this happens, the expression of a homologous gene may "mask" expression of a gene of interest, artificially skewing the results for that

35 particular gene.

Matsubara and Okubo proposed using partial cDNA sequences to establish expression profiles of genes which could be used in functional analyses of the human genome. Matsubara and Okubo warned against using random priming, as

it creates multiple unique DNA fragments from individual mRNAs and may thus skew the analysis of the number of particular mRNAs per library. They sequenced randomly selected members from a 3'-directed cDNA library and established the frequency of appearance of the various ESTs. They proposed comparing lists of ESTs from various cell types to classify genes. Genes expressed in many different cell types were labeled housekeepers and those selectively expressed in certain cells were labeled cell-specific genes, even in the absence of the full sequence of the gene or the biological activity of the gene product.

The present invention avoids the drawbacks of the prior art by providing a method to quantify the relative abundance of multiple gene transcripts in a given biological specimen by the use of high-throughput sequence-specific analysis of individual RNAs and/or their corresponding cDNAs.

The present invention offers several advantages over current protein discovery methods which attempt to isolate individual proteins based upon biological effects. The method of the instant invention provides for detailed diagnostic comparisons of cell profiles revealing numerous changes in the expression of individual transcripts.

The instant invention provides several advantages over

25 current subtraction methods including a more complex
library analysis (106 to 107 clones as compared to 103
clones) which allows identification of low abundance
messages as well as enabling the identification of messages
which either increase or decrease in abundance. These

30 large libraries are very routine to make in contrast to the
libraries of previous methods. In addition, homologues can
easily be distinguished with the method of the instant
invention.

This method is very convenient because it organizes a large quantity of data into a comprehensible, digestible format. The most significant differences are highlighted by electronic subtraction. In depth analyses are made more convenient.

The present invention provides several advantages over previous methods of electronic analysis of CDNA. The method is particularly powerful when more than 100 and preferably more than 1,000 gene transcripts are analyzed.

5 In such a case, newslow-frequency transcripts are the frequency transcripts are the frequency discovered and tissue styped.

used directly as a diagnostic profile or to identify disease-specific genesafor the development of more classic diagnostic approaches to the development of more classic This process is defined as gene transcript frequency analysis. The resulting quantitative analysis of the gene transcripts is defined as comparative gene transcript.

5 24 1 1 1 2 2 2 3 3 2 SUMMARY OF THE INVENTION OF THE 12 2

The interest of the 25 to 10 to 10 stagement of the allogical according to the

The invention is a method of analyzing a specimen containing gene transcripts comprising the steps of (a) producing a library of biological sequences; (b) generating a set of transcript sequences, where each of the transcript 20 sequences in said set is indicative of a different one of a the biological sequences of the library; (c) processing the transcript sequences in a programmed computer (in which a database of reference transcript sequences indicative of reference sequences is stored), to generate an identified sequence value for each of the transcript sequences, where each said identified sequence value is indicative of sequence annotation and a degree of match between one of the biological sequences of the library and at least one of the reference sequences; and (d) processing each said identified sequence value to generate final data values indicative of the number of times each identified sequence value is present in the library.

The invention also includes a method of comparing two specimens containing gene transcripts. The first specimen is is processed as described above. The second specimen is used to produce a second library of biological sequences, which is used to generate a second set of transcript sequences, where each of the transcript sequences in the

second set is indicative of one of the biological sequences of the second library. Then the second set of transcript sequences is processed in a programmed computer to generate a second set of identified sequence values, namely the 5 further identified sequence values, each of which is indicative of a sequence annotation and includes a degree of match between one of the biological sequences of the second library and at least one of the reference sequences. The further identified sequence values are processed to 10 generate further final data values indicative of the number of times each further identified sequence value is present in the second library. The final data values from the second first specimen and the further identified sequence values from the second specimen are processed to generate ratios at 15 of transcript sequences, which indicate the differences in the number of gene transcripts between the two specimens.

In a further embodiment, the method includes quantifying the relative abundance of mRNA in a biological specimen by (a) isolating a population of mRNA transcripts from a biological specimen; (b) identifying genes from which the mRNA was transcribed by a sequence-specific method; (c) determining the numbers of mRNA transcripts corresponding to each of the genes; and (d) using the mRNA transcript numbers to determine the relative abundance of mRNA transcripts within the population of mRNA transcripts.

Also disclosed is a method of producing a gene transcript image analysis by first obtaining a mixture of mRNA, from which cDNA copies are made. The cDNA is inserted into a suitable vector which is used to transfect 30 suitable host strain cells which are plated out and permitted to grow into clones, each cone representing a A representative population of clones unique mRNA. transfected with cDNA is isolated. Each clone in the population is identified by a sequence-specific method 35 which identifies the gene from which the unique mRNA was transcribed. The number of times each gene is identified to a clone is determined to evaluate gene transcript abundance. The genes and their abundances are listed in order of abundance to produce a gene transcript image.

In a further embodiment, the relative abundance of the gene transcripts in one cell type or tissue is compared with the relative abundance of gene transcript numbers in a second cell type or tissue in order to identify the characteristics 5 differences and similarities, we and an accommon by in a further embodiment, the method includes a system for analyzing a library of biological sequences including a means for receiving a set of transcript sequences, where each of the transcript sequences is indicative of a second 10 different one of the biological sequences of the library; and a means for processing the transcript sequences in a computer systemmin which a database of reference transcript sequences indicative of reference sequences is stored, wherein the computer is programmed with software for the second 15 generating an identified sequence value for each of the transcript sequences, where each said identified sequence value is indicative of a sequence annotation and the degree of match between a different one of the biological sequences of the library and at least one of the reference 20 sequences, and for processing each said identified sequence value to generate final data values indicative of the

In essence, the invention is a method and system for quantifying the relative abundance of gene transcripts in a biological specimen. The invention provides a method for comparing the gene transcript image from two or more different biological specimens in order to distinguish between the two specimens and identify one or more genes which are differentially expressed between the two specimens. Thus, this gene transcript image and its comparison can be used as a diagnostic. One embodiment of the method generates high-throughput sequence-specific analysis of multiple RNAs or their corresponding cDNAs: a gene transcript image. Another embodiment of the method produces the gene transcript imaging analysis by the use of high-throughput cDNA sequence analysis. In addition, two or more gene transcript images can be compared and used to detect or diagnose a particular biological state, disease,

number of times each identified sequence value is present

in the library.

I have to the a the appearance its total and the second of makes with

or condition which is correlated to the relative abundance of gene transcripts in a given cell or population of cells.

MEDICAL THE TABLES AND DRAWINGS CONTROL OF THE TABLES AND DRAWINGS

To gove were expensived. All 4.19 TABLES are position to described

5 Codes Table 1 presents a detailed explanation of the letter codes utilized in Tables 2-5.0 Codes to the letter transcripts. Table 2 lists the one hundred most common gene transcripts. This a partial list of isolates from the code HUVEC CDNA library prepared and sequenced as described

below. The left-hand column refers to the sequence's order of abundance in this table. The next column labeled to "number" is the clone number of the first HUVEC sequence identification reference matching the sequence in the read "entry" column number. Isolates that have not been

sequenced are not present in Table 2. The next column, labeled "N", indicates the total number of cDNAs which have the same degree of match with the sequence of the reference transcript in the "entry" column.

The column labeled "entry" gives the NIH GENBANK locus
name, which corresponds to the library sequence numbers.
The "s" column indicates in a few cases the species of the
reference sequence. The code for column "s" is given in
Table 1. The column labeled "descriptor" provides a plain
English explanation of the identity of the sequence
corresponding to the NIH GENBANK locus name in the "entry"
column.

Table 3 is a comparison of the top fifteen most abundant gene transcripts in normal monocytes and activated macrophage cells.

Table 4 is a detailed summary of library subtraction analysis summary comparing the THP-1 and human macrophage cDNA sequences. In Table 4, the same code as in Table 2 is used. Additional columns are for "bgfreq" (abundance number in the subtractant library), "rfend" (abundance number in the target library) and "ratio" (the target abundance number divided by the subtractant abundance number). As is clear from perusal of the table, when the abundance number in the subtractant library is "0", the

target abundance number is divided by 0.05. This is a way of obtaining a result (not possible dividing by 0) and distinguishing the result from ratios of subtractant to numbers of 1.2 180 and 10.00 and 12.00 and 12.0

Table 5 is the computer program, written in source code, for generating gene transcript subtraction profiles.

Table 6 is a partial listing of database entries used in the electronic northern blot analysis as provided by the present invention.

4.2. BRIEF DESCRIPTION OF THE DRAWINGS

Figure 1 is a chart summarizing data collected and stored regarding the library construction portion of sequence preparation and analysis.

· Production of Sept. 1994 Center of the Sept. 1994 Production of the Sept. 1995 Production of the Sep

10

25

15 <u>Figure 2</u> is a diagram representing the sequence of operations performed by "abundance sort" software in a class of preferred embodiments of the inventive method.

Figure 3 is a block diagram of a preferred embodiment of the system of the invention.

Figure 4 is a more detailed block diagram of the bioinformatics process from new sequence (that has already been sequenced but not identified) to printout of the transcript imaging analysis and the provision of database subscriptions.

5. <u>DETAILED DESCRIPTION OF THE INVENTION</u>

The present invention provides a method to compare the relative abundance of gene transcripts in different biological specimens by the use of high-throughput sequence-specific analysis of individual RNAs or their corresponding cDNAs (or alternatively, of data representing other biological sequences). This process is denoted herein as gene transcript imaging. The quantitative analysis of the relative abundance for a set of gene transcripts is denoted herein as "gene transcript image analysis" or "gene transcript frequency analysis". The present invention allows one to obtain a profile for gene transcription in any given population of cells or tissue from any type of organism. The invention can be applied to

obtain a profile of a specimen consisting of a single cell (for clones of assingle cell), or of many cells, or of tissue more complex than a single cell and containing multiple cell types, such as liver. 55 Finders The invention has significant advantages in the fields of diagnostics, toxicology and pharmacology, to name a few. A highly sophisticated diagnostic test can be performed on the ill patient in whom a diagnosis has not been made. A biological specimen consisting of the patient's fluids or ilo tissues is obtained, and the gene transcripts are isolated and expanded to the extent necessary to determine their identity. Optionally, the gene transcripts can be an about converted to cDNA. A sampling of the gene transcripts are subjected to sequence-specific analysis and quantified. 15 These gene transcript sequence abundances are compared against reference database; sequence abundances including. normal data sets for diseased and healthy patients. The an

For example, gene transcript frequency analysis can be used to differentiate normal cells or tissues from diseased cells or tissues, just as it highlights differences between normal monocytes and activated macrophages in Table 3.

patient has the disease(s) with which the patient's data

set most closely correlates.

20

In toxicology, a fundamental question is which tests
are most effective in predicting or detecting a toxic
effect. Gene transcript imaging provides highly detailed
information on the cell and tissue environment, some of
which would not be obvious in conventional, less detailed
screening methods. The gene transcript image is a more
powerful method to predict drug toxicity and efficacy.
Similar benefits accrue in the use of this tool in
pharmacology. The gene transcript image can be used
selectively to look at protein categories which are
expected to be affected, for example, enzymes which
detoxify toxins.

In an alternative embodiment, comparative gene transcript frequency analysis is used to differentiate between cancer cells which respond to anti-cancer agents and those which do not respond. Examples of anti-cancer

agents are tamoxifen, vincristine, vinblastine,
podophyllotoxins, etoposide, tenisposide, cisplatin,
biologic response modifiers such as interferon, Il-2, GMCSF, enzymes, hormones and the like. This method also
provides a means for sorting the gene transcripts by
functional category. In the case of cancer cells,
transcription factors or other essential regulatory
molecules are very important categories to analyze across
different libraries.

In yet another embodiment, comparative gene transcript frequency analysis is used to differentiate between control liver cells and liver cells isolated from patients treated with experimental drugs like FIAU to distinguish between pathology caused by the underlying disease and that caused by the drug.

In yet another embodiment, comparative gene transcript frequency analysis is used to differentiate between brain tissue from patients treated and untreated with lithium.

In a further embodiment, comparative gene transcript
20 frequency analysis is used to differentiate between
cyclosporin and FK506-treated cells and normal cells.

In a further embodiment, comparative gene transcript frequency analysis is used to differentiate between virally infected (including HIV-infected) human cells and uninfected human cells. Gene transcript frequency analysis is also used to rapidly survey gene transcripts in HIV-resistant, HIV-infected, and HIV-sensitive cells. Comparison of gene transcript abundance will indicate the success of treatment and/or new avenues to study.

In a further embodiment, comparative gene transcript frequency analysis is used to differentiate between bronchial lavage fluids from healthy and unhealthy patients with a variety of ailments.

In a further embodiment, comparative gene transcript frequency analysis is used to differentiate between cell, plant, microbial and animal mutants and wild-type species. In addition, the transcript abundance program is adapted to permit the scientist to evaluate the transcription of one gene in many different tissues. Such comparisons could

product and point mutants which produce a less abundant or otherwise different message. Such mutations can affect basic biochemical and pharmacological processes, such as mineral nutrition and metabolism, and can be isolated by means known to those skilled in the art. Thus, crops with improved yields, pest resistance and other factors can be developed.

frequency analysis is used for an interspecies comparative analysis which would allow for the selection of better pharmacologic animal models. In this embodiment, humans and other animals (such as a mouse), or their cultured cells are treated with a specific test agent. The relative sequence abundance of each cDNA population is determined. If the animal test system is a good model, homologous genes in the animal cDNA population should change expression similarly to those in human cells. If side effects are detected with the drug, a detailed transcript abundance analysis will be performed to survey gene transcript changes. Models will then be evaluated by comparing basic physiological changes.

In a further embodiment, comparative gene transcript frequency analysis is used in a clinical setting to give a 25 highly detailed gene transcript profile of a patient's cells or tissue (for example, a blood sample). In particular, gene transcript frequency analysis is used to give a high resolution gene expression profile of a diseased state or condition.

In the preferred embodiment, the method utilizes high-throughput cDNA sequencing to identify specific transcripts of interest. The generated cDNA and deduced amino acid sequences are then extensively compared with GENBANK and other sequence data banks as described below.

The method offers several advantages over current protein discovery by two-dimensional gel methods which try to identify individual proteins involved in a particular biological effect. Here, detailed comparisons of profiles of activated and inactive cells reveal numerous changes in

the expression of individual transcripts. After it is determined if the sequence is an "exact" match, similar or a non-match, the sequence is entered into a database.

Next, the numbers of copies of cDNA corresponding to each gene are tabulated. Although this can be done slowly and arduously, if at all, by human hand from a printout of all entries, a computer program is a useful and rapid way to tabulate this information. The numbers of cDNA copies (optionally divided by the total number of sequences in the data set) provides a picture of the relative abundance of transcripts for each corresponding gene. The list of represented genes can then be sorted by abundance in the cDNA population. A multitude of additional types of comparisons or dimensions are possible and are exemplified below.

An alternate method of producing asgene transcript image includes the steps of obtaining a mixture of test. mRNA and providing a representative array of unique probes whose sequences are complementary to at least some of the 20 test mRNAs. Next, a fixed amount of the test mRNA is added to the arrayed probes. The test mRNA is incubated with the probes for a sufficient time to allow hybrids of the test mRNA and probes to form. The mRNA-probe hybrids are detected and the quantity determined. The hybrids are identified by their location in the probe array. The 25 quantity of each hybrid is summed to give a population number. Each hybrid quantity is divided by the population number to provide a set of relative abundance data termed a gene transcript image analysis.

30

6. EXAMPLES

The examples below are provided to illustrate the subject invention. These examples are provided by way of illustration and are not included for the purpose of limiting the invention.

35

6.1. TISSUE SOURCES AND CELL LINES

For analysis with the computer program claimed herein, biological sequences can be obtained from virtually any

source. Most popular are tissues obtained from the human body. Tissues can be obtained from any organ of the body, any age donor, any abnormality or any immortalized cell line. Immortal cell lines may be preferred in some 5 minstances because of their purity of cell type; other tissue samples invariably include mixed cell types. A special technique is available to take a single cell (for example, a brain cell) and harness the cellular machinery to grow up sufficient cDNA for sequencing by the techniques 10 sand analysis described herein (cf. U.S. Patent Nos. 5,021,335 and 5,168,038, which are incorporated by and and , reference) ... The examples given therein utilized the trans following immortalized cell lines: monocyte-like U-937 rcells, activated macrophage-like THP-1 cells, induced to a control of the contro 15 vascular endothelial cells (HUVEC cells) and mast cell-like HMC-1 cells. The second restriction of the s The U-937 cell line is a human histiocytic lymphoma cell line with monocyte characteristics, established from malignant cells obtained from the pleural effusion of a 20 patient with diffuse histiocytic lymphoma (Sundstrom, C. and Nilsson, K. (1976) Int. J. Cancer 17:565). U-937 is one of only a few human cell lines with the morphology, cytochemistry, surface receptors and monocyte-like characteristics of histiocytic cells. These cells can be 25 induced to terminal monocytic differentiation and will express new cell surface molecules when activated with supernatants from human mixed lymphocyte cultures. Upon this type of in vitro activation, the cells undergo morphological and functional changes, including 30 augmentation of antibody-dependent cellular cytotoxicity (ADCC) against erythroid and tumor target cells (one of the principal functions of macrophages). Activation of U-937 cells with phorbol 12-myristate 13-acetate (PMA) in vitro stimulates the production of several compounds, including 35 prostaglandins, leukotrienes and platelet-activating factor (PAF), which are potent inflammatory mediators. 937 is a cell line that is well suited for the identification and isolation of gene transcripts associated with normal monocytes.

The HUVEC cell line is a normal, homogeneous, well characterized, early passage endothelial cell culture from human umbilical vein (Cell Systems Corp., 12815 NE 124th Street, Kirkland, WA 98034). Only gene transcripts from 5 induced, or treated, HUVEC cells were sequenced. One batch of 1 X 108 cells was treated for 5 hours with 1 U/ml rIL-1b and 100 ng/ml E.coli lipopolysaccharide (LPS) endotoxin prior to harvesting. A separate batch of 2 X 108 cells was treated at confluence with 4 U/ml TNF and 2 U/ml interferon-gamma (IFN-gamma) prior to harvesting. THP-1 is a human leukemic cell line with distinct monocytic characteristics. This cell line was derived from the blood of a 1-year-old boy with acute monocytic leukemia (Tsuchiya, S. et al. (1980) Int. J. Cancer: 171-76). The 15 following cytological and cytochemical criteria were used to determine the monocytic nature of the cell line: 1) the presence of alpha-naphthyl butyrate esterase activity which could be inhibited by sodium fluoride; 2) the production of lysozyme; 3) the phagocytosis of latex particles and be sensitized SRBC (sheep red blood cells); and 4) the ability of mitomycin C-treated THP-1 cells to activate Tlymphocytes following ConA (concanavalin A) treatment. Morphologically, the cytoplasm contained small azurophilic granules and the nucleus was indented and irregularly 25 shaped with deep folds. The cell line had Fc and C3b receptors, probably functioning in phagocytosis. THP-1 cells treated with the tumor promoter 12-o-tetradecanovlphorbol-13 acetate (TPA) stop proliferating and differentiate into macrophage-like cells which mimic native 30 monocyte-derived macrophages in several respects. Morphologically, as the cells change shape, the nucleus becomes more irregular and additional phagocytic vacuoles appear in the cytoplasm. The differentiated THP-1 cells also exhibit an increased adherence to tissue culture plastic.

HMC-1 cells (a human mast cell line) were established from the peripheral blood of a Mayo Clinic patient with mast cell leukemia (Leukemia Res. (1988) 12:345-55). The cultured cells looked similar to immature cloned murine

mast cells, contained histamine, and stained positively for chloroacetate esterase, amino caproate esterase, eosinophil major basic protein (MBP) and tryptase. The HMC-1 cells have, however, lost the ability to synthesize normal IgE receptors. HMC-1 cells also possess a 10;16 translocation, present in cells initially collected by leukophoresis from the patient and not an artifact of culturing. Thus, HMC-1 cells are a good model for mast cells.

CONSTRUCTION OF CONA LIBRARIES

10 For inter-library comparisons, the libraries must be prepared in similar manners. Certain parameters appear to be particularly important to control. One such parameter is the method of isolating mRNA. It is important to use the same conditions to remove DNA and heterogeneous nuclear RNA from comparison libraries. Size fractionation of cDNA must be carefully controlled. The same vector preferably should be used for preparing libraries to be compared. At the very least, the same type of vector (e.g., unidirectional vector) should be used to assure a valid comparison. A unidirectional vector may be preferred in order to more easily analyze the output.

It is preferred to prime only with oligo dT
unidirectional primer in order to obtain one only clone per
mRNA transcript when obtaining cDNAs. However, it is
recognized that employing a mixture of oligo dT and random
primers can also be advantageous because such a mixture
results in more sequence diversity when gene discovery also
is a goal. Similar effects can be obtained with DR2
(Clontech) and HXLOX (US Biochemical) and also vectors from
Invitrogen and Novagen. These vectors have two
requirements. First, there must be primer sites for
commercially available primers such as T3 or M13 reverse
primers. Second, the vector must accept inserts up to 10
kB.

It also is important that the clones be randomly sampled, and that a significant population of clones is used. Data have been generated with 5,000 clones; however, if very rare genes are to be obtained and/or their relative

abundance determined, as many as 100,000 clones from a many single library may need to be sampled. Assize fractionation of cDNA also must be carefully controlled a Alternately, plaques can be selected, rather than clones.

- 5: Besides the Uni-ZAP™ vector system by Stratagene disclosed below, it is now believed that other similarly unidirectional vectors also can be used. For example, it is believed that such vectors include but are not limited to DR2: (Clontech) and HXLOX: (U.S.: Biochemical).
- 100 Preferably, the details of library construction (as shown in Figure 1) are collected and stored in a database for later retrieval relative to the sequences being compared. Fig. 1 shows important information regarding the library collaborator or cell or cDNA supplier,
- pretreatment, biological source, culture, mRNA preparation and cDNA construction. Similarly detailed information about the other steps is beneficial in analyzing sequences and libraries in depth.
- RNA must be harvested from cells and tissue samples
 and cDNA libraries are subsequently constructed. cDNA
 libraries can be constructed according to techniques known
 in the art. (See, for example, Maniatis, T. et al. (1982)
 Molecular Cloning, Cold Spring Harbor Laboratory, New
 York). cDNA libraries may also be purchased. The U-937
 cDNA library (catalog No. 937207) was obtained from
 Stratagene, Inc., 11099 M. Torrey Pines Rd., La Jolla, CA

The THP-1 cDNA library was custom constructed by Stratagene from THP-1 cells cultured 48 hours with 100 nm TPA and 4 hours with 1 µg/ml LPS. The human mast cell HMC-1 cDNA library was also custom constructed by Stratagene from cultured HMC-1 cells. The HUVEC cDNA library was custom constructed by Stratagene from two batches of induced HUVEC cells which were separately processed.

92037.

Essentially, all the libraries were prepared in the same manner. First, poly(A+)RNA (mRNA) was purified. For the U-937 and HMC-1 RNA, cDNA synthesis was only primed with oligo dT. For the THP-1 and HUVEC RNA, cDNA synthesis was primed separately with both oligo dT and random

hexamers, and the two cDNA libraries were treated separately. Synthetic adaptor oligonucleotides were ligated onto cDNA ends enabling its insertion into the Uni-Zap^m vector system (Stratagene), allowing high efficiency of unidirectional (sense orientation) lambda library construction and the convenience of a plasmid system with blue, white color selection to detect clones with cDNA insertions. Finally, the two libraries were combined into a single library by mixing equal numbers of bacteriophage.

or antibody probes and the pBluescript® phagemid (Stratagene) can be rapidly excised in vivo. The phagemid allows the use of a plasmid system for easy insert characterization, sequencing, site-directed mutagenesis, the creation of unidirectional deletions and expression of fusion proteins. The custom-constructed library phage particles were infected into E. coli host strain XL1-Blue® (Stratagene), which has a high transformation efficiency, increasing the probability of obtaining rare, under-

6.3. ISOLATION OF CDNA CLONES

Survey of the survey great the great the

The second second second second second

The phagemid forms of individual cDNA clones were obtained by the in vivo excision process, in which the host bacterial strain was coinfected with both the lambda

25 library phage and an f1 helper phage. Proteins derived from both the library-containing phage and the helper phage nicked the lambda DNA, initiated new DNA synthesis from defined sequences on the lambda target DNA and created a smaller, single stranded circular phagemid DNA molecule

30 that included all DNA sequences of the pBluescript® plasmid and the cDNA insert. The phagemid DNA was secreted from the cells and purified, then used to re-infect fresh host cells, where the double stranded phagemid DNA was produced. Because the phagemid carries the gene for beta-lactamase,

35 the newly-transformed bacteria are selected on medium containing ampicillin.

Phagemid DNA was purified using the Magic Minipreps™
DNA Purification System (Promega catalogue #A7100. Promega

WO 95/20681 PCT/US95/01160 ·

Corp., 2800 Woods Hollow Rd., Madison, WI 53711). This small-scale process provides a simple and reliable method for lysing the bacterial cells and rapidly isolating purified phagemid DNA using a proprietary DNA-binding resin. The DNA was eluted from the purification resin already prepared for DNA sequencing and other analytical manipulations.

Phagemid DNA was also purified using the QIAwell-8
Plasmid Purification System from QIAGEN® DNA Purification
System (QIAGEN Inc., 9259 Eton Ave., Chattsworth, CA
91311). This product line provides a convenient, rapid and reliable high-throughput method for lysing the bacterial cells and isolating highly purified phagemid DNA using QIAGEN anion-exchange resin particles with EMPORE™ membrane technology from 3M in a multiwell format. The DNA was eluted from the purification resin already prepared for DNA sequencing and other analytical manipulations.

An alternate method of purifying phagemid has recently become available. It utilizes the Miniprep Kit (Catalog 20 No. 77468, available from Advanced Genetic Technologies Corp., 19212 Orbit Drive, Gaithersburg, Maryland). This kit is in the 96-well format and provides enough reagents for 960 purifications. Each kit is provided with a recommended protocol, which has been employed except for 25 the following changes. First, the 96 wells are each filled with only 1 ml of sterile terrific broth with carbenicillin at 25 mg/L and glycerol at 0.4%. After the wells are inoculated, the bacteria are cultured for 24 hours and lysed with 60 μ l of lysis buffer. A centrifugation step 30 (2900 rpm for 5 minutes) is performed before the contents of the block are added to the primary filter plate. optional step of adding isopropanol to TRIS buffer is not routinely performed. After the last step in the protocol, samples are transferred to a Beckman 96-well block for 35 storage.

Another new DNA purification system is the WIZARD product line which is available from Promega (catalog No. A7071) and may be adaptable to the 96-well format.

6.4. SEQUENCING OF CONA CLONES

Date: The cDNA inserts from random isolates of the U-937 and THP-1 libraries; were sequenced in part. Methods for DNA sequencing are well known in the art. soConventional se 5 enzymatic methods employ DNA polymerase Klenow fragment. Sequenase or Taq polymerase to extend DNA chains from an a oligonucleotide primer annealed to the DNA template of interest. Methods have been developed for the use of both single- and double-stranded templates. The chain was a great 10 termination reaction products are usually electrophoresed on urea-acrylamide gels and are detected wither by the discourse autoradiography (for radionuclide-labeled precursors) or by fluorescence (for fluorescent-labeled precursors). Recent improvements in mechanized reaction preparation, sequencing 15 and analysis using the fluorescent detection method have permitted expansion in the number of sequences that can be determined per day (such as the Applied Biosystems 373 and 377 DNA sequencer, Catalyst 800) : Currently with the system as described, read lengths range from 250 to 400 bases and are clone dependent. Read length also varies with the length of time the gel is run. In general, the shorter runs tend to truncate the sequence. A minimum of only about 25 to 50 bases is necessary to establish the identification and degree of homology of the sequence. Gene transcript imaging can be used with any sequence-2.5 specific method, including, but not limited to hybridization, mass spectroscopy, capillary electrophoresis

6.5. HOMOLOGY SEARCHING OF CDNA CLONE AND DEDUCED PROTEIN (and Subsequent Steps)

and 505 gel electrophoresis.

30

Using the nucleotide sequences derived from the cDNA clones as query sequences (sequences of a Sequence Listing), databases containing previously identified sequences are searched for areas of homology (similarity).

35 Examples of such databases include Genbank and EMBL. We next describe examples of two homology search algorithms that can be used, and then describe the subsequent computer-implemented steps to be performed in accordance with preferred embodiments of the invention.

implemented steps of the invention, the word "library" denotes a set (or population) of biological specimen nucleic acid sequences. A "library" can consist of cDNA sequences, or the like, which characterize a biological specimen. The biological specimen can consist of cells of a single human cell type (or can be any of the other above-mentioned types of specimens). We contemplate that the sequences in a library have been determined so as to accurately represent or characterize a biological specimen (for example, they can consist of representative).

specimen. (for example, they can consist of representative; cDNA sequences from clones of RNA taken from a single human cell). We have seen and the second of the second o

In the following description of the computer—
implemented steps of the invention, the expression
"database" denotes a set of stored data which represent a
collection of sequences, which in turn represent a
collection of biological reference materials. For example,
a database can consist of data representing many stored

20 cDNA sequences which are in turn representative of human
cells infected with various viruses, cells of humans of
various ages, cells from different mammalian species, and
so on.

In preferred embodiments, the invention employs a computer programmed with software (to be described) for performing the following steps:

- (a) processing data indicative of a library of cDNA sequences (generated as a result of high-throughput cDNA sequencing or other method) to determine whether each sequence in the library matches a DNA sequence of a reference database of DNA sequences (and if so, identifying the reference database entry which matches the sequence and indicating the degree of match between the reference sequence and the library sequence) and assigning an identified sequence value based on the sequence annotation and degree of match to each of the sequences in the
 - (b) for some or all entries of the database, tabulating the number of matching identified sequence

library;

values in the library (Although this can be done by human hand from a printout of all entries, we prefer to perform this step using computer software to be described below.), thereby generating a set of final data values or "abundance 5 numbers"; cand agree rule to the graphs of warding first king /s.co.(c) if the libraries are different sizes, dividing or each abundance number by the total number of sequences in the library, to obtain a relative abundance number for each identified sequence value (i.e., a relative abundance of cost 10 each gene stranscript) and Discharge and teach and the log Small The list of identified sequence values (or genes when corresponding thereto) can then be sorted by abundance in the cDNA population. Asmultitude of additional types of comparisons or dimensions (are possible. 6). Least me to the sea sec. 15 For example (to be described below in greater detail) steps (a) and (b) can be repeated for two different of the feet libraries (sometimes referred to as a "target" library and a "subtractant" library). Then, for each identified a second sequence value (or gene transcript), a "ratio" value is and 20 obtained by dividing the abundance number (for that $0 \le \varepsilon$ identified sequence value) for the target library, by the abundance number (for that identified sequence value) for the subtractant library. In fact, subtraction may be carried out on multiple 25 libraries. It is possible to add the transcripts from the or several libraries (for example, three) and then to divide them by another set of transcripts from multiple libraries (again, for example, three). Notation for this operation may be abbreviated as (A+B+C) / (D+E+F), where the capital letters each indicate an entire library. Optionally the abundance numbers of transcripts in the summed libraries may be divided by the total sample size before subtraction. 2 Unlike standard hybridization technology which permits a single subtraction of two libraries, once one has 35 processed a set or library transcript sequences and stored them in the computer, any number of subtractions can be performed on the library. For example, by this method, ratio values can be obtained by dividing relative abundance

values in a first library by corresponding values in a second library and vice versa.

In variations on step (a), the library consists of an uncleotide sequences derived from cDNA clones. Examples of databases which can be searched for areas of homology (similarity) in step (a) include the commercially available databases known as Genbank (NIH) EMBL (European Molecular Biology Labs), Germany), and GENESEQ (Intelligenetics, Mountain View, California).

implement step (a) is the algorithm which can be used to implement step (a) is the algorithm described in the paper by D.J. Lipman and W.R. Pearson, entitled "Rapid and Sensitive Protein Similarity Searches," Science, 227:1435 (1985). In this algorithm, the homologous regions are searched in a two-step manner. In the first step, the highest homologous regions are determined by calculating a matching score using a homology score table. The parameter "Ktup" is used in this step to establish the minimum window size to be shifted for comparing two sequences. Ktup also sets the number of bases that must match to extract the highest homologous region among the sequences. In this step, no insertions or deletions are applied and the homology is displayed as an initial (INIT) value.

In the second step, the homologous regions are aligned to obtain the highest matching score by inserting a gap in order to add a probable deleted portion. The matching score obtained in the first step is recalculated using the homology score Table and the insertion score Table to an optimized (OPT) value in the final output.

DNA homologies between two sequences can be examined graphically using the Harr method of constructing dot matrix homology plots (Needleman, S.B. and Wunsch, C.O., J. Mom. Biol 48:443 (1970)). This method produces a two-dimensional plot which can be useful in determining regions of homology versus regions of repetition.

However, in a class of preferred embodiments, step (a) is implemented by processing the library data in the commercially available computer program known as the INHERIT 670 Sequence Analysis System, available from

Applied Biosystems Inc. (Foster City, California), including the software known as the Factura software (also available from Applied Biosystems Inc.). The Factura ... program preprocesses each library sequence to "edit out" 5 portions thereof which are not likely to be of interest, such as the vector used to prepare the library. Additional sequences which can be edited out or masked (ignored by the search tools) include but are not limited to the polyA tail and repetitive GAG and CCC sequences. A low-end search 10 program can be written to mask out such "low-information". sequences, or programs such as BLAST can ignore the lowdinformation sequences. The week the rest time early In the algorithm implemented by the INHERIT 670 core Sequence Analysis System, the Pattern Specification 15 Language (developed by TRW Inc.) is used to determine regions of homology. "There are three parameters that determine how INHERIT analysis runs sequence comparisons: window size, window offset and error tolerance. Window size specifies the length of the segments into which the 20 query sequence is subdivided. Window offset specifies where to start the next segment [to be compared], counting from the beginning of the previous segment. Error tolerance specifies the total number of insertions, deletions and/or substitutions that are tolerated over the 25 specified word length. Error tolerance may be set to any integer between 0 and 6. The default settings are window tolerance=20, window offset=10 and error tolerance=3." INHERIT Analysis Users Manual, pp.2-15. Version 1.0, Applied Biosystems, Inc., October 1991.

Using a combination of these three parameters, a database (such as a DNA database) can be searched for sequences containing regions of homology and the appropriate sequences are scored with an initial value. Subsequently, these homologous regions are examined using dot matrix homology plots to determine regions of homology versus regions of repetition. Smith-Waterman alignments can be used to display the results of the homology search. The INHERIT software can be executed by a Sun computer system programmed with the UNIX operating system.

Search alternatives to INHERIT include the BLAST program, GCG (available from the Genetics Computer Group, WI) and the Dasher program (Temple Smith, Boston (1941) (1994) University; Boston, MA). Nucleotide sequences can be -5 searched against Genbank, EMBL or custom databases such as GENESEQ (available from Intelligenetics, Mountain View, CA) or other databases for genes. In addition, we have searched some sequences against our own in-house database. in the preferred embodiments, the transcript sequences are 10 analyzed by the INHERIT software for best conformance with a reference gene transcript to assign a sequence identifier and assigned the degree of homology, which together are the identified sequence value and are input into, and further processed by, a Macintosh personal computer (available from 15 Apple) programmed with an "abundance sort; and subtraction analysis" computer program (to be described below) ... Prior to the abundance sort and subtraction analysis program (also denoted as the "abundance sort" program), identified sequences from the cDNA clones are assigned 20 value (according to the parameters given above) by degree of match according to the following categories: "exact" matches (regions with a high degree of identity), homologous human matches (regions of high similarity, but not "exact" matches), homologous non-human matches (regions of high similarity present in species other than human), or non matches (no significant regions of homology to previously identified nucleotide sequences stored in the form of the database). Alternately, the degree of match can be a numeric value as described below. 30

With reference again to the step of identifying matches between reference sequences and database entries, protein and peptide sequences can be deduced from the nucleic acid sequences. Using the deduced polypeptide sequence, the match identification can be performed in a manner analogous to that done with cDNA sequences. A protein sequence is used as a query sequence and compared to the previously identified sequences contained in a database such as the Swiss/Prot, PIR and the NBRF Protein database to find homologous proteins. These proteins are

initially scored for homology using a homology score Table (Orcutt, B.C. and Dayoff, M.O. Scoring Matrices, PIR Report MAT - 0285 (February 1985)) resulting in an INIT score. The homologous regions are aligned to obtain the 5.5 highest, matching scores by inserting a gap which adds and approbable deleted portion. The matching score sisma on appropriate the state of th strecalculated using the homology score Table and the cinsertion score Table resulting in an optimized (OPT) score. Even in the absence of knowledge of the property of 10 greading frame of angisolated sequence, the above-described protein homology search may be performed by searching all 3 m reading frames. We great the complete protection which temp Peptide and protein sequence homologies can also be ascertained using the INHERIT 670 Sequence Analysis System 15 in an analogous way to that used in DNA sequence the sequence homologies. Pattern Specification Language and parameter windows are used to search protein databases for sequences containing regions of homology which are scored with an initial value. Subsequent display in a dot-matrix homology 20 plot shows regions of homology versus regions of repetition. Additional search tools that are available to use on pattern search databases include PLsearch Blocks (available from Henikoff & Henikoff, University of Washington, Seattle), Dasher and GCG. Pattern search 25 databases include, but are not limited to, Protein Blocks (available from Henikoff & Henikoff, University of Washington, Seattle), Brookhaven Protein (available from the Brookhaven National Laboratory, Brookhaven, MA), PROSITE (available from Amos Bairoch, University of Geneva, 30 Switzerland), ProDom (available from Temple Smith, Boston University), and PROTEIN MOTIF FINGERPRINT (available from University of Leeds, United Kingdom).

The ABI Assembler application software, part of the INHERIT DNA analysis system (available from Applied Biosystems, Inc., Foster City, CA), can be employed to create and manage sequence assembly projects by assembling data from selected sequence fragments into a larger sequence. The Assembler software combines two advanced computer technologies which maximize the ability to

assemble sequenced DNA fragments into Assemblages, a special grouping of data where the relationships between sequences are shown by graphic overlap, alignment and statistical views. The process is based on the

Meyers-Kececioglu model of fragment assembly (INHERIT Assembler User's Manual, Applied Biosystems, Inc., Foster City, CA), and uses graph theory as the foundation of a very rigorous multiple sequence alignment engine for assembling DNA sequence fragments. Other assembly programs that can be used include MEGALIGN (available from DNASTAR Inc., Madison, WI), Dasher and STADEN (available from Roger Staden, Cambridge, England).

Next, with reference to Fig. 2, we describe in more detail the "abundance sort" program which implements above—

15 mentioned "step (b) " to tabulate the number of sequences of the library which match each database entry (the "abundance number" for each database entry).

Fig. 2 is a flow chart of a preferred embodiment of the abundance sort program. A source code listing of this embodiment of the abundance sort program is set forth in Table 5. In the Table 5 implementation, the abundance sort program is written using the FoxBASE programming language commercially available from Microsoft Corporation. Although FoxBASE was the program chosen for the first iteration of this technology, it should not be considered limiting. Many other programming languages, Sybase being a particularly desirable alternative, can also be used, as will be obvious to one with ordinary skill in the art. The subroutine names specified in Fig. 2 correspond to subroutines listed in Table 5.

With reference again to Fig. 2, the "Identified Sequences" are transcript sequences representing each sequence of the library and a corresponding identification of the database entry (if any) which it matches. In other words, the "Identified Sequences" are transcript sequences representing the output of above-discussed "step (a)."

Fig. 3 is a block diagram of a system for implementing the invention. The Fig. 3 system includes library generation unit 2 which generates a library and asserts an

output stream of transcript sequences indicative of the biological sequences comprising the library. Programmed processor 4 receives the data stream output from unit 2 and processes this data in accordance with above-discussed.

5 "step (a) " to generate the Identified Sequences. Processor 4 can be a processor programmed with the commercially available computer program known as the INHERIT 670 Sequence Analysis System and the commercially available computer program known as the Factura program (both available from Applied Biosystems Inc.) and with the UNIX operating system.

Still with reference to Fig. 3, the Identified Sequences are loaded into processor 6 which is programmed with the abundance sort program. Processor 6 generates the Final Transcript sequences indicated in both Figs. 2 and 3. Fig. 4 shows a more detailed block diagram of a planned relational computer system, including various searching techniques which can be implemented, along with an assortment of databases to query against.

With reference to Fig. 2, the abundance sort program first performs an operation known as "Tempnum" on the Identified Sequences, to discard all of the Identified Sequences except those which match database entries of selected types. For example, the Tempnum process can select Identified Sequences which represent matches of the following types with database entries (see above for definition): "exact" matches, human "homologous" matches, "other species" matches representing genes present in species other than human), "no" matches (no significant regions of homology with database entries representing previously identified nucleotide sequences), "I" matches (Incyte for not previously known DNA sequences), or "X" matches (matches ESTs in reference database). eliminates the U, S, M, V, A, R and D sequence (see Table 1 for definitions).

The identified sequence values selected during the "Tempnum" process then undergo a further selection (weeding out) operation known as "Tempred." This operation can, for

example, discard all identified sequence values representing matches with selected database entries. The identified sequence values selected during the "Tempred" process are then classified according to library, 5 during the "Tempdesig" operation. It is contemplated that the "Identified Sequences" can represent sequences from a single library, or from two or more libraries. Consider first the case that the identified sequence values represent sequences from a single library. In this 10 case, all the identified sequence values determined during "Tempred" undergo sorting in the "Templib" operation, further sorting in the "Libsort" operation, and finally additional sorting in the "Temptarsort" operation. For her example, these three sorting operations can sort the Bo identified sequences in order of decreasing "abundance and number" (to generate a list of decreasing abundance numbers, each abundance number corresponding to a unique identified sequence entry, or several lists of decreasing abundance numbers, with the abundance numbers in each list 20 corresponding to database entries of a selected type) with redundancies eliminated from each sorted list. In this case, the operation identified as "Cruncher" can be bypassed, so that the "Final Data" values are the organized transcript sequences produced during the "Temptarsort" 25 operation.

We next consider the case that the transcript sequences produced during the "Tempred" operation represent sequences from two libraries (which we will denote the "target" library and the "subtractant" library). For example, the target library may consist of cDNA sequences from clones of a diseased cell, while the subtractant library may consist of cDNA sequences from clones of the diseased cell after treatment by exposure to a drug. For another example, the target library may consist of cDNA sequences from clones of a cell type from a young human, while the subtractant library may consist of cDNA sequences from clones of the same cell type from the same human at different ages.

In this case, the "Tempdesig" operation routes all transcript sequences representing the target library for processing in accordance with "Templib" (and then "Libsort" and "Temptarsort"), and routes all transcript sequences 5 representing the subtractant library for processing in accordance with "Tempsub" (and then "Subsort" and "Tempsubsort"). For example, the consecutive "Templib," "Libsort," and "Temptarsort" sorting operations sort identified sequences from the target library in order of decreasing abundance number (to generate a list of decreasing abundance numbers, each abundance number corresponding to a database entry, or several lists of decreasing abundance numbers, with the abundance numbers in each list corresponding to database entries of a selected 15 type) with redundancies eliminated from each sorted list. The consecutive "Tempsub," "Subsort," and "Tempsubsort" sorting operations sort identified sequences from the subtractant library in order of decreasing abundance number (to generate a list of decreasing abundance numbers, each abundance number corresponding to a database entry, or several lists of decreasing abundance numbers, with the abundance numbers in each list corresponding to database entries of a selected type) with redundancies eliminated from each sorted list.

25 The transcript sequences output from the "Temptarsort" operation typically represent sorted lists from which a histogram could be generated in which position along one (e.g., horizontal) axis indicates abundance number (of target library sequences), and position along another 30 (e.g., vertical) axis indicates identified sequence value (e.g., human or non-human gene type). Similarly, the transcript sequences output from the "Tempsubsort" operation typically represent sorted lists from which a histogram could be generated in which position along one 35 (e.g., horizontal) axis indicates abundance number (of subtractant library sequences), and position along another (e.g., vertical) axis indicates identified sequence value (e.g., human or non-human gene type).

The transcript sequences (sorted lists) output from the Tempsubsort and Temptarsort sorting operations are combined during the operation identified as "Cruncher." The "Cruncher" process identifies pairs of corresponding 5 target and subtractant abundance numbers (both representing the same identified sequence value), and divides one by the other to generate a "ratio" value for each pair of corresponding abundance numbers, and then sorts the ratio values in order of decreasing ratio value. The data output 10 from the "Cruncher" operation (the Final Transcript sequence in Fig. 2) is typically a sorted list from which a histogram could be generated in which position along one: axis indicates the size of a ratio of abundance numbers (for corresponding identified sequence values from target 15 and subtractant libraries) and position along another axis indicates identified sequence value (e.g., gene type).

Preferably, prior to obtaining a ratio between the two library abundance values, the Cruncher operation also divides each ratio value by the total number of sequences in one or both of the target and subtractant libraries. The resulting lists of "relative" ratio values generated by the Cruncher operation are useful for many medical, scientific, and industrial applications. Also preferably, the output of the Cruncher operation is a set of lists, each list representing a sequence of decreasing ratio values for a different selected subset (e.g. protein family) of database entries.

In one example, the abundance sort program of the invention tabulates for a library the numbers of mRNA

30 transcripts corresponding to each gene identified in a database. These numbers are divided by the total number of clones sampled. The results of the division reflect the relative abundance of the mRNA transcripts in the cell type or tissue from which they were obtained. Obtaining this final data set is referred to herein as "gene transcript image analysis." The resulting subtracted data show exactly what proteins and genes are upregulated and downregulated in highly detailed complexity.

HOW THE RESERVENCE OF THE PROPERTY OF THE PROP

Table 2 is an abundance table listing the various gene transcripts in an induced HUVEC library. The transcripts are listed in order of decreasing abundance. This computerized sorting simplifies analysis of the tissue and speeds identification of significant new proteins which are specific to this cell type. This type of endothelial cell lines tissues of the cardiovascular system, and the more that is known about its composition, particularly in

response to activation, the more choices of protein targets become available to affect in treating disorders of this tissue, such as the highly prevalent atherosclerosis.

6.7.: (MONOCYTE-CELL AND MAST-CELL CONA-LIBRARIES)

Tables 3 and 4 show truncated comparisons of two

libraries. In Tables 3 and 4 the "normal monocytes" are
the HMC-1 cells, and the "activated macrophages" are the
THP-1 cells pretreated with PMA and activated with LPS.

Table 3 lists in descending order of abundance the most
abundant gene transcripts for both cell types. With only

20 15 gene transcripts from each cell type, this table permits quick, qualitative comparison of the most common transcripts. This abundance sort, with its convenient side-by-side display, provides an immediately useful research tool. In this example, this research tool

discloses that 1) only one of the top 15 activated
macrophage transcripts is found in the top 15 normal
monocyte gene transcripts (poly A binding protein); and 2)
a new gene transcript (previously unreported in other
databases) is relatively highly represented in activated

macrophages but is not similarly prominent in normal macrophages. Such a research tool provides researchers with a short-cut to new proteins, such as receptors, cell-surface and intracellular signalling molecules, which can serve as drug targets in commercial drug screening

programs. Such a tool could save considerable time over that consumed by a hit and miss discovery program aimed at identifying important proteins in and around cells, because those proteins carrying out everyday cellular functions and

represented as steady state mRNA are quickly eliminated from further characterization.

This illustrates how the gene transcript profiles change with altered cellular function. Those skilled in the standard the blochemical composition of cells also changes with other functional changes such as cancer, including cancer's various stages, and exposure to toxicity. A gene transcript subtraction profile such as in Table 3 is useful as a first screening tool for such gene expression and protein studies.

6.8. SUBTRACTION ANALYSIS OF NORMAL MONOCYTE-CELL AND ACTIVATED MONOCYTE CELL CONA LIBRARIES

Once the cDNA data are in the computer, the computer program as disclosed in Table 5 was used to obtain ratios; 15 of all the gene transcripts in the two libraries discussed in Example 6.7, and the gene transcripts were sorted by the descending values of their ratios. If a gene transcript is not represented in one library, that gene transcript's abundance is unknown but appears to be less than 1. As an approximation -- and to obtain a ratio, which would not be possible if the unrepresented gene were given an abundance of zero -- genes which are represented in only one of the two libraries are assigned an abundance of 1/2. Using 1/2 for unrepresented clones increases the relative importance 25 of "turned-on" and "turned-off" genes, whose products would be drug candidates. The resulting print-out is called a. subtraction table and is an extremely valuable screening method, as is shown by the following data.

Table 4 is a subtraction table, in which the normal monocyte library was electronically "subtracted" from the activated macrophage library. This table highlights most effectively the changes in abundance of the gene transcripts by activation of macrophages. Even among the first 20 gene transcripts listed, there are several unknown gene transcripts. Thus, electronic subtraction is a useful tool with which to assist researchers in identifying much more quickly the basic biochemical changes between two cell types. Such a tool can save universities and pharmaceutical companies which spend billions of dollars on

research valuable time and laboratory resources at the early discovery stage and can speed up the drug development cycle; , which: in .turn permits researchers to set up drug est screening programs much earlier. Thus, this research tool 5 provides a way to get new drugs to the public faster and more deconomically. We shall be a like to a protection of the state of these Also, such a subtraction table can be obtained for n patient diagnosis. An individual patient sample (such as monocytes obtained from a biopsy or blood sample) can be 10 compared with data provided herein to diagnose conditions associated with macrophage activation. Table 4 uncovered many new gene transcripts (labeled Incyte clones). Note that many genes are turned on in the activated macrophage (i.e., the monocyte had a 0 in the 15 bgfreq column). This screening method is superior to other screening techniques, such as the western blot, which are incapable of uncovering such a multitude of discrete new gene transcripts.

The subtraction-screening technique has also uncovered 20 a high number of cancer gene transcripts (oncogenes rho, ETS2, rab-2 ras, YPT1-related, and acute myeloid leukemia mRNA) in the activated macrophage. These transcripts may be attributed to the use of immortalized cell lines and are inherently interesting for that reason. This screening 25 technique offers a detailed picture of upregulated transcripts including oncogenes, which helps explain why anti-cancer drugs interfere with the patient's immunity mediated by activated macrophages. Armed with knowledge gained from this screening method, those skilled in the art can set up more targeted, more effective drug screening programs to identify drugs which are differentially effective against 1) both relevant cancers and activated macrophage conditions with the same gene transcript profile; 2) cancer alone; and 3) activated macrophage 35 conditions.

Smooth muscle senescent protein (22 kd) was upregulated in the activated macrophage, which indicates that it is a candidate to block in controlling inflammation.

6.9. SUBTRACTION ANALYSIS OF NORMAL LIVER CELLS AND HEPATITIS INFECTED LIVER CELL CDNA LIBRARIES

In this example, rats are exposed to hepatitis virus and maintained in the colony until they show definite signs of hepatitis. Of the rats diagnosed with hepatitis, one half of the rats are treated with a new anti-hepatitis agent (AHA). Liver samples are obtained from all rats before exposure to the hepatitis virus and at the end of AHA treatment or no treatment. In addition, liver samples can be obtained from rats with hepatitis just prior to AHA treatment.

The liver tissue is treated as described in Examples 6.2 and 6.3 to obtain mRNA and subsequently to sequence cDNA. The cDNA from each sample are processed and analyzed for abundance according to the computer program in Table 5. The resulting gene transcript images of the cDNA provide detailed pictures of the baseline (control) for each animal and of the infected and/or treated state of the animals. cDNA data for a group of samples can be combined into a group summary gene transcript profile for all control samples, all samples from infected rats and all samples from AHA-treated rats.

Subtractions are performed between appropriate individual libraries and the grouped libraries. For 25 individual animals, control and post-study samples can be subtracted. Also, if samples are obtained before and after AHA treatment, that data from individual animals and treatment groups can be subtracted. In addition, the data for all control samples can be pooled and averaged. The 25 control average can be subtracted from averages of both post-study AHA and post-study non-AHA cDNA samples. If pre- and post-treatment samples are available, pre- and post-treatment samples can be compared individually (or electronically averaged) and subtracted.

These subtraction tables are used in two general ways. First, the differences are analyzed for gene transcripts which are associated with continuing hepatic deterioration or healing. The subtraction tables are tools to isolate the effects of the drug treatment from the underlying basic pathology of hepatitis. Because hepatitis affects many

parameters, additional liver toxicity has been difficult to detect with only blood tests for the usual enzymes. The gene transcript profile and subtraction provides a much more complex biochemical picture which researchers have the subtraction analyze such difficult problems.

Second, the subtraction tables provide a tool for identifying clinical markers, individual proteins or other biochemical determinants which are used to predict and/or evaluate a clinical endpoint, such as disease, improvement 10 due to the drug, and even additional pathology due to the drug. The subtraction tables specifically highlight genes which are turned on or off. Thus, the subtraction tables provide a first screen for a set of gene transcript candidates for use as clinical markers. Subsequently, 15 electronic subtractions of additional cell and tissue libraries reveal which of the potential markers are in fact found in different cell and tissue libraries. Candidate gene transcripts found in additional libraries are removed from the set of potential clinical markers. Then, tests of 20 blood or other relevant samples which are known to lack and have the relevant condition are compared to validate the selection of the clinical marker. In this method, the particular physiologic function of the protein transcript need not be determined to qualify the gene transcript as a 25 clinical marker.

6.10. ELECTRONIC NORTHERN BLOT

One limitation of electronic subtraction is that it is difficult to compare more than a pair of images at once. Once particular individual gene products are identified as relevant to further study (via electronic subtraction or other methods), it is useful to study the expression of single genes in a multitude of different tissues. In the lab, the technique of "Northern" blot hybridization is used for this purpose. In this technique, a single cDNA, or a probe corresponding thereto, is labeled and then hybridized against a blot containing RNA samples prepared from a multitude of tissues or cell types. Upon autoradiography,

the pattern of expression of that particular gene, one at a time, can be quantitated in all the included samples.

In contrast, a further embodiment of this invention is the computerized form of this process, termed here

5 "electronic northern blot." In this variation, a single gene is queried for expression against a multitude of prepared and sequenced libraries present within the database. In this way, the pattern of expression of any single candidate gene can be examined instantaneously and effortlessly. More candidate genes can thus be scanned, leading to more frequent and fruitfully relevant discoveries. The computer program included as Table 5 includes a program for performing this function, and Table 6 is a partial listing of entries of the database used in the electronic northern blot analysis.

6.11. PHASE I CLINICAL TRIALS

the transfer of the

Based on the establishment of safety and effectiveness in the above animal tests, Phase I clinical tests are undertaken. Normal patients are subjected to the usual 20 preliminary clinical laboratory tests. In addition, appropriate specimens are taken and subjected to gene transcript analysis. Additional patient specimens are taken at predetermined intervals during the test. specimens are subjected to gene transcript analysis as described above. In addition, the gene transcript changes noted in the earlier rat toxicity study are carefully evaluated as clinical markers in the followed patients. Changes in the gene transcript analyses are evaluated as indicators of toxicity by correlation with clinical signs 30 and symptoms and other laboratory results. In addition, subtraction is performed on individual patient specimens and on averaged patient specimens. The subtraction analysis highlights any toxicological changes in the treated patients. This is a highly refined determinant of 35 toxicity. The subtraction method also annotates clinical markers. Further subgroups can be analyzed by subtraction analysis, including, for example, 1) segregation by

occurrence and type of adverse effect; and 2) segregation by dosage.

6.12. GENE TRANSCRIPT IMAGING ANALYSIS IN CLINICAL STUDIES

A gene transcript imaging analysis (or multiple gene

5 transcript imaging analyses) is a useful tool in other
clinical studies. For example, the differences in gene
transcript imaging analyses before and after treatment can
be assessed for patients on placebo and drug treatment.
This method also effectively screens for clinical markers

to follow in clinical use of the drug.

6.13. COMPARATIVE GENE TRANSCRIPT ANALYSIS BETWEEN SPECIES

The subtraction method can be used to screen cDNA libraries from diverse sources. For example, the same cell types from different species can be compared by gene transcript analysis to screen for specific differences, such as in detoxification enzyme systems. Such testing aids in the selection and validation of an animal model for the commercial purpose of drug screening or toxicological testing of drugs intended for human or animal use. When the comparison between animals of different species is shown in columns for each species, we refer to this as an interspecies comparison, or zoo blot.

Embodiments of this invention may employ databases such as those written using the FoxBASE programming
25 language commercially available from Microsoft Corporation. Other embodiments of the invention employ other databases, such as a random peptide database, a polymer database, a synthetic oligomer database, or a oligonucleotide database of the type described in U.S. Patent 5,270,170, issued
30 December 14, 1993 to Cull, et al., PCT International Application Publication No. WO 9322684, published November 11, 1993, PCT International Application Publication No. WO 9306121, published April 1, 1993, or PCT International Application Publication No. WO 9119818, published December 26, 1991. These four references (whose text is incorporated herein by reference) include teaching which

may be applied in implementing such other embodiments of the present invention. Past 1

All references referred to in the preceding text are hereby expressly incorporated by reference herein.

Swi5 of a Various modifications and variations of the described method and system of the invention will be apparent to those skilled in the art without departing from the scope and spirit of the invention. Although the invention has been described in connection with specific preferred 10 embodiments, it should be understood that the invention as 10-6 claimed should not be unduly limited to such specific embodiments. Online 1 June 10 June $C(\cdot)$

Hartelon Pickler Gorden

Wilder Green Committee Committe

to refer to authors when the stop

1. 47

17 36 36 1 36 36 760

5 18.000 5 18.000

į

5

.

7

÷ .

...

्रा १५५ १८ च्या

15 000

f = 0 is

 $\mathcal{L}_{i} = \mathcal{O}_{i} \mathcal{A}_{i}$ 1. e(E 15 000

30 000 10 115 11 650

3 10%

··....

41.4

1000

; 27. 18.

 $\mathbf{t}_{i}^{-1}(\xi_{i}) = \mathbf{t}_{i}^{-1}(\xi_{i})$

1.

150

. If the . .

TABLE 1	Localization (Z)	fic N = Nuclear ue specific C = Cytopli K = Cytosk E = Cell su Z = Intracy M = Mitoch S = Secrety U = Unknow X = Other X = Other 1 = Do pri 2 = Primar 3 = Full l 3 = Full l 5 = Tissue 5 = Obtain	43, 45 100 100 100 100 100 100 100 10
	Distribution (F)	C = Non-specific V = Cell/tissue specific V = Unknown K = Ape P = Pig D = Dog V = Bovine B = Rabit R = Rat R = Rat R = Rat C = Chicken F = Amphibian I = Invertebrate Z = Protozoan G = Fungi	
	Designations (D)	E = Exact H = Homologous O = Other species N = No match D = Noncoding gene U = Nonreadable R = Repetitive DNA A = Poly-A only V = Vector only V = Wetcor only V = Mitochondrial DNA S = Skip I = Match Incyte clone X = EST match Library (L) U = U937 M = HMC T = THP-1 H = HUVEC S = Spleen L = Lung Y = T & B cell A = Adenoid	

TABLE 2

Clone numbers 15000 through 20000 Libraries: HUVEC Arranged by ABUNDANCE
Total clones analyzed: 5000
319 genes, for a total of 1713 Clones

		Number N				
	number	N c	entry	8	descriptor	
_	1.4.	 Mail Chester of All Control All Control 	<u>-</u>			
1	15365	67	HSRPL41		Riboptn L41	
2	₅ .15004	65	NCY015004		INCYTE 015004	
3	15638	63 ₃₅ .55 (· NCY015638		INCYTE 015638	
4	15390	50	NCY015390		INCYTE 015390	
5	C25193	47	HSFIB1		Fibronectin	
6	%15220	47	RRRPL9	R	Ribonto I.9	
7	15280	<i>3</i> ′ - 47	NCY015280	••	דארעייי הובספה	
8	15583	33	M62060		FOT BUCKOO (TOD)	
9	15662	31	HSACTCCP		Action commo	
10	15026	20	NCVO15026		ACCIN, Gamma	
11	15279	24	NCIOIDO20		INCITE UI5U26	
12	15027	1. 1,.24	NOVOLEGOS		EII 1-alpha	
13	15027	23	NC1015027		INCYTE 015027	
14	15100	20	NC1015033		INCYTE 015033	
16	15150	20	NC1012138		INCYTE 015198	
15	12003	20	HSCOLLI		Collagenase	
10	15261	19	NCY015221		INCYTE 015221	
1/	15263	19	NCY015263		INCYTE 015263	
10	15290 ,	19	NCY015290	· · · · · · · · · · · · · · · · · · ·	INCYTE 015290	
73	15350	18	NCY015350	5.0	INCYTE 015350	
20	15030	17	NCY015030		INCYTE 015030	
21	15234	17	NCY015234		INCYTE 015234	
22	15459	16	NCY015459		INCYTE 015459	
23	15353	15	NCY015353	• • • •	INCYTE 015353	
24	15378	15	S76965	•	Ptn kinase inhib	
25	15255	14	HUMTHYB4		Thymosin beta-4	
26	15401	14	HSLIPCR		Lipocortin I	
27	15425	14	HSPOLYAB	· .	Poly-A bp	
28	18212 -	- 14	HUMTHYMA:		Thymosin, alpha	
29	18216	14	HSMRP1	· · · · · · · · · · · · · · · · · · ·	Motility relat ptn: MRP-1:CD-9	
30	15189	13	HS18D		Interferon induc ntn 1-80	
31	15031	12	HUMFKBP		FK506 bp	
32	15306	12	HSH2AZ		Histone H2A	
33	15621	12	HUMLEC		Lectin, B-galbo, 14kha	
34	· 15789	11	NCY015789		INCYTE 015789	
35	÷ 16578	11	HSRPS11		Ribonth S11	
36	16632	11	M61984		EST HECA13 /TCD)	
37	18314	11	NCY018314		INCYTE 018314	
38	15367	10	NCY015367		INCYTE 016314	
39	15415	10	HSIFNINI		interferon indua mays	
40	15633	10	HSLDHAR		Lactate debudrossass	
41	15813	10	CHKNMHCB		C Musein beam shair n	
42	18210	10	NCVOIRZIO		Thouse Ologio	
43	18233	10	HCDDTT140		DNA malamana ==	
44	18996	10	NCVOIROR		ANA polymerase II	
45	15088	9	HIMERDI		Enmitte Disysb	
46	15714	9	NOVOIEZIA		rerritin, light chain	
47	15720	9	NOVO15700		INCITE U15714	
7 / 10	15063	7	MCIOT2/50		INCYTE 015720	
40	15101	7	NCIOT2863		INCYTE 015863	
47	10757	9	HSET		Endothelin	
50	18252	9	NCY018252 HUMALBP		INCYTE 018252	
51	15351	8	HUMALBP		Lipid bp, adipocyte	
52	15370	8	NCY015370		INCYTE 015863 Endothelin INCYTE 018252 Lipid bp, adipocyte INCYTE 015370	

TABLE 2 Con't

	Fig. 19		-		
	Section	and the second	TAB.	LE 2 Con't	
		1 1	· •		
	number	te se verific	•		
	301 2 2 2 2 2 2 2 2	💢 5 g 😘		_	A = = = 1 b
	number	N E	entry	, 5	descriptor
E 2		`	BTCIASHI	v	NADIT
53	15670 15795 16245	8	NCY015795		NADH-ubiq oxidoreductase
55	16946	. o 1432	MOVOLEDAE		INCYTE 015795 INCYTE 016245
56	18262	e e	NCTO10243		INCITE 010245
57	10202	8	, NCIU10202	17. 8 4 18 18 18 18	INCYTE 018262 Riboptn L17
58	; 18321 - 15126	7	XLRPL1BRF		
59	· 15133	ر بن 			Riboptn Ll Actin, beta
		7	NCY015245		INCYTE 015245
61	and the second s		NCY015288		INCYTE 015245
62	15294	7 ()	Hecypub NCIOI3200	LANGE THE STA	C-3-DD
63	15442	7 7 7 7	HUMLAMB		Laminin receptor, 54kDa
64	15485	ż	HSNGMRNA		Uracil DNA glycosylase
65	16646	7 7 7	NCY016646		INCYTE 016646
66			BUMBATA		Plsmnogen activ gene
67	15032 15267	6" " 1 1 1 1	HUMUB		Uracil DNA glycosylase INCYTE 016646 Plsmnogen activ gene Ubiquitin Riboptn S8
68	15267	6	HSRPS8		Riboptn S8
69	15295	6	NCY015295		INCYTE 015295
70	15458	6	RNRPS10R		Riboptn S10
71		6	RSGALEM		UDP-galactose epimerase
72	15928	6	HUMAPOJ		Apolipoptn J
73	16598	6	HUMTBBM40		Tubulin, beta
74	18218	6	NCY018218		INCYTE 018218 Hydrophobic ptn p27 INCYTE 018963 INCYTE 018997 Galactosidase A, alpha
75	18499	6	HSP27		Hydrophobic ptn p27
76		6	NCY018963		INCYTE 018963
77	18997	6	NCY018997		INCYTE 018997
78		5	HSAGALAR		Galactosidase A, alpha
79		5	NCY015475		INCITE 0154/5
80	15721	5, 5	NCY015721		INCYTE 015721
81	15865		NCY015865		INCYTE 015865
82	16270		NCY016270		INCYTE 016270
83 84	16886		NCY016886		INCYTE 016886
85	18500		NCY018500		INCYTE 018500
86	18503 19672	5 5	NCY018503		INCYTE 018503
87	15086	4	RRRPL34 XLRPL1AR	R	Riboptn L34
	15113	4	HUMIFNWRS	·F	Riboptn Lla
89	15242	4	NCY015242		tRNA synthetase, trp INCYTE 015242
90	15249	4	NCY015242		· - · - · -
91		-	NCY015377		======
92	15407		NCY015407		INCYTE 015377 INCYTE 015407
93	15473	4	NCY015473		INCYTE 015407
94	15588		HSRPS12	•	
95	15684	4	HSEF1G		Riboptn S12 Elf 1-gamma
96	15782	-	NCY015782		INCYTE 015782
97	15916		HSRPS18		Riboptn S18
98	15930	4	NCY015930		INCYTE 015930
		-	NCY016108		INCYTE 016108
100		4	NCY016133		INCYTE 016133
		÷ ,			010133

ORMAL MONONCYTE VS. **IIVATED MACROPHAGI**

Top 15 Most Abundant Genes

NORMAL

Ribosomal protein S8 homolog Ribosomal phosphoprotein Elongation factor-l alpha Ribosomal protein L7 Ferritin H chain **Beta-Globin**

Ribosomal protein S20 homolog **Nucleoplasmin**

ransferrin receptor

Poly-A binding protein

ranslationally controlled tumor ptn Ribosomal protein S25

Signal recognition particle SRP9 Histone HZĂ.Z

Ribosomal protein Ke-3

Macrophage inflammatory protein-Interleukin-I beta nterleukin-8

Lymphocyte activation gene Elongation factor-l alpha

Beta actin

Rantes T-cell specific protein Poly A binding protein

Tumor Necrosis Factor-alpha Osteopontin; nephropontin NCYTE clone 01 1050

Adenylate cyclase (yeast homolog) Cu/Zn superoxide dismutase

NGF-related B cell activation molecule Protease Nexin-I, glial-derived

TABLE 4

: 4 .

Libraries: THP-1
Subtracting: HMC
Sorted by ABUNDANCE
Total clones analyzed: 7375

1057 genes, for a total of 2151 clones

Part 2 hall Company				
number (Mentry 4	(v. 18) descriptor 22	bgfreq	rfend	ratio
10022 HUMIL1	"The IL"1-beta - 15th 1241 to 18	0	131	262.00
10036 HSMDNCF	F IL-8	0	119	238.00
10089 HSLAG10	CDN Lymphocyte activ gene	Ŏ	71	142.00
10060 HUMTCSM		Ö	23	46.000
10003 HUMMIP1		3	121	40.333
10689 HSOP	Osteopontin	ŏ	20	40.000
11050 NCY0110		ŏ	17	34.000
10937 HSTNFR	TNF-alpha	ŏ	17	34.000
10176 HSSOD		ŏ	14	28.000
10886 HSCDW40	• · · · · · · · · · · · · · · · · · · ·	ŏ	10	20.000
10186 HUMAPR	Early resp PMA-induc	ŏ	9	18.000
10967 HUMGDN	PN-1, glial-deriv	Ö	9	18.000
11353 NCY0113		ŏ	8	16.000
10298 NCY0102	298 INCYTE 010298	Ö	7	14.000
10215 HUM4COL	LA Collagenase, type IV	Ö .	6	12.000
10276 NCY0102		ŏ	6	12.000
10488 NCY0104		Ö	6	12.000
11138 NCY0111	138 INCYTE 011138	Ō	6	12.000
10037 : HUMCAPP	PRO Adenylate cyclase	ī	10	10.000
10840 HUMADCY	Adenylate cyclase	Ó	5	10.000
10672 HSCD44E	Cell adhesion glptn	0	5	10.000
12837 HUMCYCL		0	5	10.000
10001 NCY0100	001 INCYTE 010001	Ó	5	10.000
10005 NCY0100		0	5	10.000
10294 NCY0102		0	5	10.000
10297 NCY0102		0	5	10.000
10403 NCY0104		0	5	10.000
10699 NCY0106		0	5	10.000
10966 NCY0109		0	5	10.000
12092 NCY0120		0	5	10.000
12549 HSRHOB	Oncogene rho	0	5	10.000
10691 HUMARF1		0	4	8.000
12106 HSADSS	Adenylosuccinate syntheta	se O	4	8.000
10194 HSCATHL		0	4	8.000
10479 CLMCYCA		0	4	8.000
10031 NCY0100		0	4	8.000
10203 NCY01020		0	4 .	8.000
10288 NCY01028	88 INCYTE 010288	0	4	8.000
10372 NCY0103		0	4	8.000
10471 - NCY0104: 10484 NCY01048		0	4	8.000
		0	4	8.000
10859 NCY01089		Ō	4	8.000
		Ō	4	8.000
	•	0	4	8.000
		0	4	8.000
		0	4	8.000
10133 HSI1RAP 10516 HUMP2A	IL-1 antagonist	0	4	8.000
	Phosphatase, regul 2A	0	4	8.000
	TNF-induc response NA HB15 gene; new Ig	0	4	8.000
		0	3	6.000
10788 NCY00173 10033 NCY01003		0	3	6.000
10035 NCY01003		0	3	6.000
10035 NCY01003		0	3	6.000
10236 NCY01023	_ :	0	3	6.000
10383 NCY01023		0	3	6.000
TOTOS NCIUIUSE	83 INCYTE 010383	0	3	6.000

TABLE 4 Con't

18 - 12 30 3 5 180 - 180 - 180 - 180 - 180 - 180 - 180 - 180 - 180 - 180 - 180 - 180 - 180 - 180 - 180 - 180 - 180 - 180 - 180

			TABLE 4 C	on't			
138 L WA	1012 121 - 16 - 31	Marie .					
2. F. St. 18	ভর্কার্যানীর কর্মা সংক্রান্ত করি বিশ্বাস্থা						
Grant, and Arthritis							
number	entry s	descri	ptor		bgfreq	rfend	ratio
10450 TOWN	NCY010450	TNOVTE	010450		0	2	6 000
	NCY010470		010470		Ö	3 3	6.000 6.000
10504	NCY010504	INCYTE			Ö	3	6.000
	NCY010507			3640 P 85 F			
10598	NCY010598	TNCVTE	010598	44 Pen 244 C			6.000
	NCY010779	INCYTE			0	3	6.000
10909	NCY010909	INCYTE			ŏ	3	6.000
10976	NCY010976	INCYTE			Ö	3	6.000
10985	NCY010976 NCY010985	INCYTE			Ö	3	6.000
11052	NCY011052	INCYTE			ŏ	3	6.000
11068	NCY011068	INCYTE			ŏ	3	6.000
11134	NCY011134	INCYTE			ŏ	333333333333333333333333333333333333333	6.000
	NCY011136	INCYTE	011136		Ö	3	6.000
11191 5500	NCY011191	INCYTE			ŏ	3	6.000
11219	'NCY011219	INCYTE	011219		ŏ	3	6.000
11386	NCY011386	INCYTE	011386		ō	3	6.000
11403	NCY011403	INCYTE	011403		O	3	6.000
11460	NCY011460	INCYTE	011460		0	3	6.000
11618	NCY011618	INCYTE	011618		0	3 3 3 3	6.000
11686	NCY011686	INCYTE	011686		0	3	6.000
12021	NCY012021	INCYTE	012021		0	3 3 3	6.000
	NCY012025	INCYTE	012025		0	3	6.000
12320	NCY012320	INCYTE	012320		0	3	6.000
12330	NCY012330 ·	INCYTE	012330		0	3 3	6.000
12853	NCY012853	INCYTE	012853		0	3	6.000
14386	NCY014386	INCYTE	014386		0	3	6.000
14391	NCY014391	INCYTE	014391		0	3	6.000

TABLE 5

```
The man of the second second of the second o
                              1 3 to 35 8 4
          * Master menu for BUBTRACTION output
         SET TALK OFF
                                                                                                                                                                                                                                                                                           ेक्किया ३८८ 🕝 ह
          SET EXACT ON A CONTRACT OF SERVICE
          SET TYPENHEAD TO OGRAFIANT STATE OF STATE AS A STATE OF
                                                                                                                                                                                                                                                                                                                                                                                                                  . : ::
                                                                                                                                                                                                                                                                                                                                      10 mg/ g/ sp
                                                                                                                                                                                                                                                                                                                                                                                                                                                      ATHER CONTRACTOR
          CLEAR .
         SET DEVICE TO SCREEN
          USE . * SmartGuy: FoxBASE+/Mac: fox files: Clones.dbf*
         STORE NUMBER TO INITIATE
         GO BOTTOM
         STORE NUMBER TO TERMINATE
        STORE: TO Object1
STORE: TO Object1
STORE: TO Object2
STORE: TO Object3
         STORE O TO ANAL
        STORE 0 TO HMATCH
         STORE D TO CHATCH
        STORE 0.TO IMPORTAGE CONTROL OF THE STORY OF
         STORE 1 TO BAIL
                                                                                                                                                                                                      48 182 8732 1 4 5
       DO WHILE .T.
         Program. 1 Subtraction 2.fmt
                                                                                                                                                                                             LIFE YOU THE SEASON DESCRIPTION OF THE OWNERS OF COMMITTEE
    * Program: 1 30/11/94

* Date...: 10/11/94

* Version: FoxEASE+/Mac, revision 1.10

* Notes...: Format file Subtraction 2
SCREEN 1 TYPE 0 HEADING "Screen 1" AT 40,2 SIZE 286,492 PIXELS PONT "Geneva",9 COLOR 0,0,0,0 FIXELS 75,120 TO 178,241 STYLE 3871 COLOR 0,0,-1,24610,-1,8947

6 FIXELS 27,134 SAY "Subtraction Memu" STYLE 65536 FONT "Geneva",274 COLOR 0,0,-1,-1,-1,-1 PIXELS 117,126 GET EMATCH STYLE 65536 FONT "Chicago";12 PICTURE "6"C Exact " SIZE 15,62 CO 6 FIXELS 135,126 GET EMATCH STYLE 65536 FONT "Chicago",12 PICTURE "6"C Homologous" SIZE 15,1 FIXELS 153,126 GET CMATCH STYLE 65536 FONT "Chicago",12 PICTURE "6"C Other spc" SIZE 15,84 FIXELS 90,152 SAY "Matches:" STYLE 65536 FONT "Chicago",12 PICTURE "6"C Other spc" SIZE 15,84 FIXELS 90,152 SAY "Matches:" STYLE 65536 FONT "Chicago",12 PICTURE "6"C Incyte" SIZE 15,65 CO 6 PIXELS 252,137 GET Initiate STYLE 0 FONT "Geneva",12 SIZE 15,70 COLOR 0,0,-1,-1,-1,-1 PIXELS 252,236 GET terminate STYLE 0 FONT "Geneva",12 SIZE 15,70 COLOR 0,0,-1,-1,-1,-1 FIXELS 252,235 SAY "Include clones" STYLE 3536 FONT "Geneva",12 SIZE 15,70 COLOR 0,0,-1,-1,-1,-1,-1 FIXELS 252,215 SAY ">> STYLE 65536 FONT "Geneva",14 COLOR 0,0,-1,-1,-1,-1,-1 FIXELS 198,126 GET PIT STYLE 35536 FONT "Geneva",14 COLOR 0,0,-1,-1,-1,-1,-1 FIXELS 198,126 GET PIT STYLE 3871 COLOR 0,0,-1,-25600,-1,-1

6 PIXELS 90,9 TO 181,109 STYLE 3871 COLOR 0,0,-1,-25600,-1,-1

6 PIXELS 81,296 SAY "Background: STYLE 65536 FONT "Chicago",12 PICTURE "6"R Overall; Function" STZE 4 FIXELS 81,26 SAY "Background: STYLE 65536 FONT "Geneva",270 COLOR 0,0,-1,-1,-1,-1,-1

7 PIXELS 81,296 SAY "Background: STYLE 65536 FONT "Geneva",270 COLOR 0,0,-1,-1,-1,-1,-1

8 PIXELS 81,26 SAY "Target: STYLE 65536 FONT "Geneva",270 COLOR 0,0,-1,-1,-1,-1,-1

9 PIXELS 81,26 SAY "Target: STYLE 65536 FONT "Geneva",270 COLOR 0,0,-1,-1,-1,-1,-1

9 PIXELS 81,26 SAY "Target: STYLE 0 FONT "Geneva",270 COLOR 0,0,-1,-1,-1,-1,-1,-1

9 PIXELS 81,270 SAY "Target: STYLE 0 FONT "Geneva",270 COLOR 0,0,-1,-1,-1,-1,-1,-1
 0 FIXELS 81,26 EAY "Target:" STYLE 65536 FONT "Geneva",270 COLOR 0,0,-1,-1,-1,-1
9 FIXELS 108,20 GET target1 STYLE 0 FONT "Geneva",9 SIZE 12,79 COLOR 0,0,-1,-1,-1,-1
9 FIXELS 135,20 GET target2 STYLE 0 FONT "Geneva",9 SIZE 12,79 COLOR 0,0,-1,-1,-1,-1
10 FIXELS 162,20 GET target2 STYLE 0 FONT "Geneva",9 SIZE 12,79 COLOR 0,0,-1,-1,-1,-1
10 FIXELS 108,299 GET object1 STYLE 0 FONT "Geneva",9 SIZE 12,79 COLOR 0,0,-1,-1,-1,-1
10 FIXELS 135,299 GET object2 STYLE 0 FONT "Geneva",9 SIZE 12,79 COLOR 0,0,-1,-1,-1,-1
10 FIXELS 162,299 GET object3 STYLE 0 FONT "Geneva",9 SIZE 12,79 COLOR 0,0,-1,-1,-1,-1
10 FIXELS 276,324 GET Bail STYLE 65536 FONT "Chicago",12 FICTURE "6"R Rum,Bail out" SIZE 4112
      * EOF: Subtraction. 2.fmt
     READ
               IF Bailm2
               CLEAR
CLOSE DATABASES
                USE "SmartGuy: FoxBASE+/Mac: fox files: clones.dbf"
               SET SAFETY ON
                 SCREEN 1 OFF
               RETURN
```

```
The other was an experience of the second property of the second of the 
                    EMITTE
          STORE VAL(SYS(2)) TO STARTIME
          STORE UPPER (Target1) TO Target1 TO Target2
          STORE UPPER (Target3) TO Target3
                                                                                                                                                                                                                                                          The state of the s
           STORE UPPER (Object1) TO Object1
         STORE UPPER (Object2). TO Object2
STORE UPPER (Object3). TO Object3
STORE UPPER (Object3). TO Object3
STORE UPPER (Object3). TO Object3
         GAP # TERMINATE-INITIATE+1 PARTY OF THE WAR ARE A PARTY OF THE PARTY OF A PARTY OF A PARTY OF THE PARTY OF TH
         GO INITIATE
         COPY NEXT GAP FIELDS NUMBER, Library, D. P. Z. R. ENTRY, S. DESCRIPTOR, START, RFEND, I TO TEMPNUM
        USE TEMPNUM
        COUNT TO TOT.
        COFY TO TEMPRED FOR D='E'.OR.D='O'.OR.D='H'.OR.D='N'.OR.D='I'
      IF Ematched AND, Bhatched AND, Onatched AND, Dearched
      COPY TO TEMPDESIG
ELSE
COPY STRUCTURE TO TEMPDESIG
                                                                                                                                                                                                                                                                                                                                                                  AND THE STATE OF SEASON SERVICES OF COURSE OF THE
      USE TEMPORATE

TE PRATCHAL

LANGUAGE

LANGUAGE
                      APPEND FROM TEMPNUM FOR D='B'
                   ENDIF THE TOTAL THE STATE OF TH
                  APPEND FROM TEMPNUM FOR D='H'
                   IF Onatch=1
APPEND FROM TEMPNUM FOR D=101
                   IF Imatchel
                   APPEND FROM TEMPNIM FOR DE'I'.OR.DE'X' - SANCE AND THE SEE FOR TELLINGISC AND CORRECT
     ENDIF
   COUNT TO STARTOT A STARTOT AND THE START OF 
    COPY STRUCTURE TO TEMPLIE
  USE TEMPLIB . .
                 APPEND FROM TEMPDESIG FOR library=UPPER(target1)
                 IP target2co
               APPEND FROM TEMPDESIG FOR library=UPPER(target2)
ENDIF
IF target3<>'
APPEND FROM TEMPDESIG FOR library=UPPER(target3)
ENDIF
ENDIF
COUNT TO ANALMOT
USE TEMPDZSIG
                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                           Ten Toler Oliver to
                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                     The second of th
                                                                                                                                                                                                                                                          COPY STRUCTURE TO TEMPSUB
  USE TEMPSUB
                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                      7 . W. W. 3
               APPEND FROM TEMPDESIG FOR library=UPPER(Object1)
               IP target2co'
               APPEND FROM TEMPDESIG FOR library=UFFER(Object2)
             ENDIF
               IF target3co'
            APPEND FROM TEMPDESIG FOR library=UPPER(Object3)
             ENDIP
 COUNT TO SUBTRACTOT
 SET TALK OFF
 * COMPRESSION SUBROUTINE A ? 'COMPRESSING QUERY LIBRARY'
 USE TEMPLIE
```

the other Commences States of the money of the all the complete with the

```
emission from a received the court as as as a size and a size of the many the separate size of the second
                 SORT ON ENTRY, NUMBER TO LIBSORT, A. P. C. . TOWN D. A. M. C. GRADO C. DA SOLVETTER AND MILES.
             USE LIESORT MADDALETA COUNTY OF COUNTY OF THE STATE OF TH
             MARKI - 1
              5W2=0 1.19
                                                                                                                                           Market topics of the II of the South State of the World Continues and
             DO WHILE SW2=0 ROLL
                A IP MARKI >= IDGRNB A CONTROL OF AN ENGINEER SERVED WAS SERVED ON A SERVED OF A CONTROL OF AN A SERVED OF A CONTROL OF AN A CONTROL OF AN A CONTROL OF A CONTROL
                 ASSOCIATION OF THE CONTRACT OF THE PROPERTY OF THE STATE OF THE CONTRACT OF THE STATE OF THE STA
                      ENDIF
                                                                                                                                                                                                                                                 海绵等领品建筑的 建筑 电电流设计
            GO MARKI
            DUP = 1
          STORE ENTRY TO TESTA, THE STANDARD OF THE SHEET WAS A TO SHEET WITH THE STANDARD STA
            STORE D TO DESIGA
         SW = 0. The second of the seco
                                                                                       199
          STORE ENTRY TO TESTE AND TO AN ADDRESS OF AN ADDRESS OF AN ADDRESS OF A DESCRIPTION OF AN ADDRESS OF A DESCRIPTION OF A DESCR
                 IF TESTA - TESTB.AND.DESIGA-DESIGB
DELSTE
DUP = DUP-1
                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                              THE REPORT OF THE PARTY.
                      DUP = DUP+1
                        LOOP
                                                                                                                                          where the moving of the day of the contract
                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                  · •
                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                               I To too how they
                  MDIF
        GO MARKI
        REPLACE REEND WITH DUP
        MARK1 = MARK1+DUP '
        SW=1
LCOP
      LOOP
ENUDO. TEST
LOOP
                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                           Sec. 30 1.25 1.25 1.33.25 1.
                                                                                                                                                                                                                                                                                 en der Mille State (1994), der Mille
Mille State (1994), der State (1994),
                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                  Committee Commit
        SORT ON REPEND/D, NUMBER TO TEMPTARSORT.
      USE TEMPTARSORT
        *REPLACE ALL START WITH RPEND/IDGENE*10000
     COUNT TO TEMPTARCO
     * COMPRÉSSION SUBROUTINE B
     ? 'COMPRESSING TARGET LIBRARY'
    USE TEMPSUB
     SORT ON ENTRY, NUMBER TO SUBSORT
     use Subsort
  COUNT TO SUBGENE
REPLACE ALL REEND WITH 1
  MARK1 = 1
     6W2=0
DO WHILE SW2=0 ROLL
                  IF MARKL >= SUEGENE
                  PACK .
                 COUNT TO BUNIQUE
                  SW2=1
                LOOP .
                 ENDIF
 GO MARKI .
 DUP = 1
  STORE ENTRY TO TESTA
  STORE D TO DESIGN
  5W = 0
 DO WHILE SW-0 TEST
 BKIP
 STORE ENTRY TO TESTE
 STORE D TO DESIGN
                IF TESTA = TESTB.AND.DESIGA=DESIGB
```

```
DUF = DUF+1
                                                                            1. 18 Car 18 19
     LOOP
ENDIP
     REPLACE REEND WITH DUP
     MARKI - MARKI+DUP
     MAKALA SHE
                                                                                                                                        \{(i,1,\dots)_m\}
     ENDOO TEST TO 1 S
     ENDDO ROLL
     FORT ON REFEND/D, NUMBER TO TEMPSUBSORT
   USE TEMPSUBSORT
   *REPLACE ALL START WITH RPEND/INCEME*10000
COUNT TO TEMPSORCO COUNT TO
 FUBICAL ROUTINE
? SUBTRACTION
COPY STRUCTURE TO CRUNCHER
SELECT 2
USE TEMPSUBSORT
USE CRUNCHER
APPEND FROM TEMPTARSORT
COUNT TO BAILOUT
MARK = 0
  .T. EIIIW CO
  SELECT 1
  MARK = MARK+1
       IF MARK-BAILOUT
        EXIT
       ENDIP
 GO MARK
  STORE ENTRY TO SCANNER
 SHEECT 2
LOCATE FOR ENTRY-SCANNER
  IF FOUND() .
 STORE REEND TO BITT
 STORE REEND TO BIT?
 else ·
 STORE 1/2 TO BITT
 STORE 0 TO BIT2
BNDIF
 ENDIF
SELECT 1
 REPLACE BGFREO WITH BITT2
 REPLACE ACTUAL WITH BITT
 LOOP
 ENDO
                                      2.00
 SELECT 1
REPLACE ALL RATIO WITH RFEND/ACTUAL
7 'DOING FINAL SORT BY RATIO'
SORT ON RATIO/D, EGFREQ/D, DESCRIPTOR TO FINAL
USE PINAL
met talk off
DO CASE
CASE PTF=0
SET DEVICE TO PRINT
SET PRINT ON
EJECT ...
CASE PIF=1
SET ALTERNATE TO "Adenoid Patent Figures: Subtraction.txt"
```

```
SET ALTERNATE ON
ENDCASE
    STORE VALISYS(2)) TO FINTIME
IF FINTIME-STARTIME
STORE FINTIME+86400 TO FINTIME
ENDIF
      STORE FINTIME - STARTIME TO COMPSEC
     STORE COMPSEC/60 TO COMPMIN
    , which are as the constraint of the street field of the 1950 days.
     SET MARGIN TO 10
     61,1 SAY 'Library Subtraction Analysis' STYLE 65536 FORT 'Geneva',274 COLOR 0,0,0,-1,-1,-1
    ?? TDE()
? 'Clone numbers'
    ? Target1
    IP Target2<>'
                                                           Street of the Market of the 
Market Blood for Market of the 
Market of the Market of the Alles of the 
Market of the Market of the Market of the Market of the Market of the Market of the Market of the Market of the Market of the Market of the Market of the Market of the Market of the Market of the Market of the Market of the Market of the Market of the Market of the Market of the Market of the Market of the Market of the Market of the Market of the Market of the Market of the Market of the Market of the Market of the Market of the Market of the Market of the Market of the Market of the Market of the Market of the Market of the Market of the Market of the Market of the Market of the Market of the Market of the Market of the Market of the Market of the Market of the Market of the Market of the Market of the Market of the Market of the Market of the Market of the Market of the Market of the Market of the Market of the Market of the Market of the Market of the Market of the Market of the Market of the Market of the Market of the Market of the Market of the Market of the Market of the Market of the Market of the Market of the Market of the Market of the Market of the Market of the Market of the Market of the Market of the Market of the Market of the Market of the Market of the Market of the Market of the Market of the Market of the Market of the Market of the Market of the Market of the Market of the Market of the Market of the Market of the Market of the Market of the Market of the Market of the Market of the Market of the Market of the Market of the Market of the Market of the Market of the Market of the Market of the Market of the Market of the Market of the Market of the Market of the Market of the Market of the Market of the Market of the Market of the Market of the Market of the Market of the Market of the Market of the Market of the Market of the Market of the Market of the Market of the Market of the Market of the Market of the Market of the Market of the Market of the Market of the Market of the Market of t
    77 Target2
    PADIF
    IP Target3<>'.
    33 14, 1,
    ?? Target3
    ENDIF
    ? 'Subtracting:
    ? Object1
    IF-Object2<>'
    ?? . .
    ?? Object2
   ENDIF
    IF Object3<>'
   ?? Object3
   ENDIF
  ? 'Designations: '
  IF Smatch=0 .AND. Hmatch=0 .AND. Omatch=0 .AND. IMATCH=0 ?? 'All'
  ENDIF
IF Ematch=1
   ?? 'Exact,'
  ENDIF
  IF Hmatch=1
?? 'Human,'
  EVIDIF
  IF Ometch=1
   ?? 'Other sp.'
  ENDIF
 IF Imatch=1
??'INCYTE'
  ENDIF
· IF ANAL=1
  7 'Sorted by ABUNDANCE'
RNDIF.
IF ANAL=2
? 'Arranged by FUNCTION'
```

. ...

C .. .

Ţ.,

```
? 'Total clones représented: '
               ?? STR (TOT, 5, 0)
                ? 'Total clones analyzed: '
               ?? STR (STARTOT, 5, 0)
               ? 'Total computation time:
              ?? STR (COMPMIN, 5, 2)
               ?? | mimites!
     70'd = debignation 60'f = distribution s = location r = function s = species
(30)** ******************
         BCREEN 1 TYPE 0 HEADING "Screen 1" AT 40,2 SIZE 286,492 PIXELS FONT "Geneva",9 COLOR 0.0.0.
 DO CASE .
 CASE ANAL-1
### ANALIST MANALIST AND PROPERTY OF THE PROPE
             SCREEN 1 TYPE 0 HEADING "Screen 1" AT 40,2 SIZE 286,492 PIXELS FONT "Geneva",7 COLOR 0,0,0, list OFF fields number, D, F, Z, R, ENTRY, S, DESCRIPTOR, EGFREQ, RFEND, RATIO, I
1711
              SET PRINT OFF
              CLOSE DATABASES
       USE 'SmartGuy: FoxBASE+/Mac: fox files: clones.dbf'
( ) H ( )
             CASE. ANAL=2
 arrange/function
SET PRINT ON
             SET HEADING ON
             SCREEN 1 TYPE O HEADING 'Screen 1' AT 40,2 SIZE 285,492 PIXELS FONT 'Helvetica',268 COLOR O
                                                                                BINDING PROTEINS
             SCREEN 1 TYPE 0 READING "Screen 1: AT 40,2 SIZE 285,492 PIXELS FORT "Helvetica",265 COLOR 0
             ? 'Surface molecules and receptors:'

SCREEN 1 TYPE 0 HEADING "Screen 1" AT 40,2 SIZE 285,492 PIXELS FORT "Geneva",7 COLOR 0,0,0,
List CFF fields number, D, F; Z, R, ENTRY, S, LESCRIPTOR, BGFRED, RREND, RATIO, I FOR R='B'
             SCREWN 1 TYPE 0 HEADING "Screen 1" AT 40.2 SIZE 286.492 PIXELS .FONT "Helvetica", 265 COLOR 0
   ? 'Calcium-binding proteins:'
            SCREEN 1 TYPE 0 HEADING "Screen 1" AT 40,2 SIZE 286,492 PIXELS FONT "Geneva",7 COLOR 0,0,0, list OFF fields number, D.F.Z.R. ENTRY, S. DESCRIPTOR, BGFRED, RFEND, RATIO, I FOR RE"C"
            SCREEN 1 TYPE 0 HEADING "Screen 1" AT 40,2 SIZE 286,492 PIXELS PONT "Helvetica",265 COLOR 0
            ? 'Ligands and effectors:'

SCREEN 1 TYPE 0 HEADING "Screen 1" AT 40,2 SIZE 286,492 PIXELS FONT "Geneva",7 COLOR 0,0,0,
List OFF fields number, D, F, Z, R, EMTRY, S, DESCRIPTOR, BGFRED, RFEND, RATIO, I FOR R='S'
            SCREEN 1 TYPE 0 HEADING 'Screen 1' AT 40,2 SIZE 286,492 FIXELS FORT "Helvetica",265 COLOR 0
             ? 'Other binding proteins:'
            SCREEN 1 TYPE 0 HEADING 'Screen 1' AT 40,2 SIZE 286,492 PIXELS FONT "Geneva",7 COLOR 0,0,0, list OFF fields number, D, F, Z, R, ENTRY, S, DESCRIPTOR, EGFRSQ, RFEND, RATIO, I FOR R='I'
            SCREEN 1 TYPE 0 HEADING 'Screen 1' AT 40,2 SIZE 286,492 PIXELS FORT 'Helvetica',268 COLOR 0
                                                                                   ONCOGENES!
             SCREEN 1 TYPE 0 HEADING 'Screen 1' AT 40,2 SIZE 286,492 PIXELS FONT 'Helvetica',265 COLOR 0
            ? 'General oncogenes:',
SCREEN 1 TYPE 0 HEADING 'Screen 1' AT.40,2 SIZE 286,492 PIXELS FONT 'Geneva',7 COLOR 0,0,0,
list OFF fields number, D, F, Z, R, ENTRY, S, DESCRIPTOR, EGFRED, RATIO, I FOR Re'O'
            SCREEN 1 TYPE 0 HEADING 'Screen 1' AT 40,2 SIZE 286,492 PIXELS FORT 'Helvetica',265 COLOR 0
            ? 'GTP-binding proteins!'
            SCREEN 1 TYPE 0 HEADING 'Screen 1' AT 40,2 SIZE 286,492 PIXELS FORT 'Geneva',7 COLOR 0,0,0,
            list OFF fields number, D, F, Z, R, EMTRY, S, DESCRIPTOR, BGFREQ, RFEND, RATIO, I FOR R='Q'
```

```
SCREEN 1 TYPE 0 HEADING "Screen 1" AT 40,2 SIZE 286,492 PIXELS FONT "Helvetica",265 COLOR 0
        ? 'Viral elements: 'Screen 1' AT 40,2 SIZE 286,492 PIXELS FONT "GENEVA", 7 COLOR 0,0,0, 1ist OFF fields number, D, F, Z, R, ENTRY, S, DESCRIPTOR, EGFREQ, RFEND, RATIO, I FOR RE'V'
        SCREEN 1 TYPE 0 HEADING *Screen 1* AT 40,2 SIZE 286,492 PIXELS FONT *Helvetica*,265 COLOR 0
         ? 'Kinases and Phosphatases:' .
        SCREEN 1 TYPE 0 HEADING "Screen 1" AT 40,2 SIZE 286,492 PIXELS FORT "Geneva",7 COLOR 0,0,0,
        list OFF fields number, D.F. Z. R. ENTRY, S. DESCRIPTOR, BOFREQ, RFEND, RATIO, I FOR Ra'Y'
        SCREEN 1 TYPE 0 HEADING "Screen 1" AT 40,2 SIZE 286,492 PIXELS PONT "Helvetica",265 COLOR 0
        ? 'Tumor-related antigens:'
        SCREEN 1 TYPE 0 HEADING "Screen 1" AT 40.2 SIZE 286,492 PIXELS FONT "Geneva",7 COLOR 0,0,0,
 list OFF fields number, D, F, Z, R, ENTRY, S, DESCRIPTOR, ECFREQ, RFEND, RATIO, I FOR R='A'
SCREEN 1 TYPE O HEADING "Screen 1" AT 40.2 SIZE 286,492 PIXELS FONT "Helvetica", 268 COLOR O
                        PROTEIN SYMPHETIC MACHINERY PROTEINS
SCREEN 1 TYPE 0 HEADING "Screen 1" AT 40,2 SIZE 286,492 PIXELS FONT "Helvetica",265 COLOR 0 Transcription and Nucleic Acid-binding proteins:

SCREEN 1 TYPE 0 HEADING "Screen 1" AT 40,2 SIZE 286,492 PIXELS FONT "Geneva",7 COLOR 0,0,0, list OFF fields number, D, F, Z, R, ENTRY, S, DESCRIPTOR, BOFREQ, RFEND, RATIO, I FOR R='D'
SCREEN 1 TYPE 0 HEADING "Screen 1" AT 40,2 SIZE 286,492 PIXELS FONT "Helvetica",265 COLOR 0
     SCREEN 1 TYPE 0 HEADING "Screen 1" AT 40.2 SIZE 286,492 PIXELS FONT "Geneva",7 COLOR 0.0.0. list OFF fields number, D. F. Z. R. ENTRY, S. DESCRIPTOR, EGFREQ, RFEND, RATIO, I FOR R='T'
     SCREEN 1 TYPE 0 HEADING "Screen 1" AT 40,2 SIZE 286,492 PIXELS FONT "Helvetica",265 COLOR 0
       ? 'Riboscmal proteins:'
SCREEN 1 TYPE 0 HEADING "Screen 1" AT 40,2 SIZE 286,492 PIXELS FONT "Geneva",7 COLOR 0,0,0,
  list OFF fields number, D; F, Z, R, ENTRY, S, DESCRIPTOR, BGFREQ, RFEND, RATIO, I FOR Re'R'
   SCREEN 1 TYPE 0 HEADING "Screen 1" AT 40.2 SIZE 286,492 PIXELS FONT "Helvetica",265 COLOR 0
       ? 'Protein processing:
       SCREEN 1 TYPE 0 HEADING "Screen 1" AT 40,2 SIZE 286,492 PIXELS FONT "Geneva",7 COLOR 0,0,0,
       list OFF fields number, D.F.Z.R. ENTRY, S. DESCRIPTOR, BGFREQ, RFEND, RATIO, I FOR Ra'L'
      SCREEN 1 TYPE 0 HEADING "Screen 1" AT 40,2 SIZE 285,492 PIXELS, FONT "Helvetica", 268 COLOR 0
                                                ENZYMES!
       SCREEN 1 TYPE 0 HEADING 'Screen 1' AT 40,2 SIZE 286,492 PIXELS FONT "Helvetica",265 COLOR 0
       ? 'Ferroproteins:'
      SCREEN 1 TYPE 0 HEADING "Screen 1" AT 40,2 SIZE 286,492 PIXELS FONT "Geneva",7 COLOR 0.0.0. list OFF fields number, D, F, Z, R, ENTRY, S, DESCRIPTOR, BGFREQ, RFEND, RATIO, I FOR R='F'
      SCREEN 1 TYPE 0 HEADING *Screen 1* AT 40,2 SIZE 285,492 FIXELS FONT "Helvetica",265 COLOR 0
      ? 'Proteases and inhibitors:'
      SCREEN 1 TYPE 0 HEADING "Screen 1" AT 40.2 SIZE 286.492 PIXELS FORT "Geneva".7 COLOR 0.0.0, list OFF fields number; D.F.Z.R.EWIRY, S.DESCRIPTOR, BGFREQ, RFEND, RATIO, I FOR Re'P'
      SCREEN 1 TYPE 0 HEADING "Screen 1" AT 40,2 SIZE 285,492 PIXELS FONT "Helvetica",265 COLOR 0
      ? 'Oxidative phosphorylation:' ... SCREEN 1 TYPE O HEADING "Screen 1" AT 40.2 SIZE 286.492 FIXELS FONT "Geneva", 7 COLOR 0,0,0,
      list OFF fields number, D, F, Z, R, ENTRY, S, DESCRIPTOR, BGFREQ, RFEND, RATIO, I FOR R='Z'
      SCREEN 1 TYPE 0 HEADING "Screen 1" AT 40,2 SIZE 286,492 PIXELS FONT "Helvetica",255 COLOR 0
      ? 'Sugar metabolism:'
      SCREEN 1 TYPE O HEADING "Screen 1" AT 40,2 SIZE 286,492 PIXELS FORT "Geneva".7 COLOR 0,0,0, list OFF fields number, D, F, Z, R, EMTRY, S, DESCRIPTOR, EGFREQ, RFEND, RATIO, I FOR R='Q'
      SCREEN 1 TYPE 0 HEADING "Screen 1" AT 40,2 SIZE 286,492 FIXELS FORT "Helvetica",265 COLOR 0
      7 'Amino acid metabolism:'
      SCREEN 1 TYPE 0 HEADING "Screen 1" AT 40,2 SIZE 286,492 PIXELS FONT "Geneva",7 COLOR 0,0,0,
```

111.7

list OFF fields number, D, F, Z, R, EMTRY, S, DESCRIPTOR, BOFREQ, RFEND, RATIO, I FOR R='M' SCREEN 1 TYPE 0. HEADING "Screen 1" AT 40,2 SIZE 286,492 PINCLS PONT "Relvelica", 253 COLOR O ? 'Nucleic acid metabolism; + SCREEN 1.TYPE 0 HEADING "Screen 1" AT 40.2 SIZE 286,492 FIXELS FONT "Geneva",7 COLOR 0,0,0, list OFF fields number, D, F, Z, R, ENTRY, S, DESCRIPTOR, BOFFEQ, RFEND, RATID, I FOR RE'N' ECREEN 1 TYPE 0 HEADING "Screen 1" AT 40,2 SIZE 285,492 PIXELS FONT "Helvetica", 265 COLOR 0 Tall 1-7 !Lipid metabolism: SECREEN 1 TYPE: 0 HEADING "Screen 1" AT 40,2 SIZE 286,492 PIXELS FONT "Geneve",7 COLOR 0,0,0, GEO. A. LIST OFF fields number, D, F, Z, R, ENTRY, S, DESCRIPTOR, BGFREQ, RFEND, RATIO, I FOR RE'W' SCREEN 1 TYPE 0 HEADING "Screen 1" AT 40,2 SIZE 286,492 PIXELS FONT "Helvetica",265 COLOR 0 COLOR 0 SCREEN 1 TYPE 0 HEADING "Screen 1" AT 40,2 SIZE 286,492 PIXELS FONT "Geneva",7 COLOR 0,0,0,0 list OFF fields number, D. F. Z. R. ENTRY, S. DESCRIPTOR, EGFREQ, RFEND, RATIO, I FOR R='E' SCREEN 1 TYPE 0 HEADING "Screen 1" AT 40,2 SIZE 285,492 PIXELS FONT "Helvetice",268 COLOR 0 The result of the second of th SCREEN 1 Type 0 HEADING "Screen 1" AT 40,2 SIZE 286,492 PIXELS FONT "Helvetica", 265 COLOR 0 SCREEN 1 TYPE 0 HEADING 'Screen 1" AT 40,2 SIZE 286,492 PIXELS FONT 'Geneva',7 COLOR 0,0,0,0 list OFF fields number, D, F, Z, R, ENTRY, S, DESCRIPTOR, BGFREQ, RFEND, RATIO, I FOR R='H' SCREEN 1 TYPE 0 HEADING "Screen 1" AT 40,2 SIZE 286,492 PIXELS FONT "Helvetica", 265 COLOR'O ? 'Structural:' SCREEN 1 TYPE 0 HEADING "Screen 1" AT 40,2 SIZE 286,492 PIXELS FONT "Geneva",7 COLOR 0,0,0, list OFF fields number, D.F. Z. R. ENTRY, 8, DESCRIPTOR, BGFRED, RETIO, I FOR R='K' SCREEN 1 TYPE 0 HEADING "Screen 1" AT 40,2 SIZE 286,492 PIXELS FONT "Helvetica",265 COLOR 0 ? 'Other clones:' SCREEN 1 TYPE 0 HEADING "Screen 1" AT 40.2 SIZE 286,492 PIXELS FONT "Geneva", 7 COLOR 0.0.0. list OFF fields number, D. F. Z. R. ENTRY, S. DESCRIPTOR, BCFREQ, RFEND, RATIO, I FOR R='X'

SCREEN 1 TYPE 0 HEADING "Screen 1" AT 40,2 SIZE 286,492 PIXELS FONT "Helvetica",265 COLOR 0 ? 'Clones of unknown function:' SCREEN 1 TYPE 0 HEADING "Screen 1" AT 40,2 SIZE 286,492 PIXELS FONT "Geneva",7 COLOR 0,0,0, list OFF fields number, D, F, Z, R, ENTRY, S, DESCRIPTOR, EGFREQ, RFEND, RATIO, I FOR Re-U-ENDCASE

DO "Test print.prg" SET PRINT OF SET DEVICE TO SCREEN CLOSE DATABASES ERASE TEMPLIB.DEF ERASE TEMPNUM.DBF ERASE TEMPDESIG.DBF SET MARGIN TO 0 CLEAR LCOP ENDDO

7.17

 $::_{\mathbb{R}^{n}} \to$

```
*Northern (single), version 11-25-94
                                           close databases
                                          SET TALK OFF
                                          SET PRINT OFF
                                         SET EXACT OFF
                                         CLEAR .
                                         STORE .
                                                                                                                     ' TO Eabject
 STORE ' TO

STORE 0 TO Numb

STORE 0 TO Zog

STORE 1 TO Bail

DO WHILE T.

Program.: Northern (single) fmt

Data...: 8/8/94

Version.: FoxBASE+/Mac, revision 1.10

Notes....: Format file Northern (single)
                                        STORE
                                                                                                                                                                                                                           ' TO Dobject
SCREEN 1 TYPE 0 HEADING "Screen 1" AT 40.2 SIZE 286,492 PIXELS PONT "Geneva",12 COLOR 0,0.0 & PIXELS 15.81 TO 46,397 STYLE 28447 COLOR 0,0.-1,-25600,-1,-1 & PIXELS 89,79 TO 192,422 STYLE 28447 COLOR 0,0.0.-25600,-1,-1 & PIXELS 115.98 SAY "Entry $!" STYLE 65536 FONT "Geneva",12 COLOR 0,0.0,-1,-1,-1 & PIXELS 115.173 GET Dobject STYLE 0 FONT "Geneva",12 SIZE 15,142 COLOR 0,0.0,-1,-1,-1 & PIXELS 145.89 SAY "Description" STYLE 65536 FONT "Geneva",12 COLOR 0,0.0,-1,-1,-1 & PIXELS 145.173 GET Dobject STYLE 0 FONT "Geneva",12 SIZE 15.241 COLOR 0,0.0,-1,-1,-1 & PIXELS 35.89 SAY "Single Northern search screen" STYLE 65536 FONT "Geneva";274 COLOR 0,0,-1,-1,-1 & PIXELS 220,162 GET Bail STYLE 65536 FONT "Chicago",12 PICTURE "&*R Continue;Bail cut" SIZE PIXELS 175.98 SAY "Clone $:" STYLE 65536 FONT "Geneva",12 COLOR 0,0,-1,-1,-1 & PIXELS 175.173 GET Numb STYLE 0 FONT "Geneva",12 SIZE 15.70 COLOR 0,0,0,-1,-1,-1 & PIXELS 175.173 GET Numb STYLE 0 FONT "Geneva",12 SIZE 65536 FONT "Geneva",12 COLOR 0,0,0,-1,-1,-1,-1 & PIXELS 80,152 SAY "Enter any ONE of the following: STYLE 65536 FONT "Geneva",12 COLOR -1,
               *** ************* (single) . fmt
   READ IF Bail=2
                           CLEAR
 screen 1 off
RETURN
ENDIF
                                USE "SmartGuy: FoxEASE+/Mac:Fox files:Lookup.dbf"
                            SET TALK ON
                                IF Exbjects'
                                                                                                                                                                                                                 Company of the
                                STORE UPPER (Bobject) to Bobject
                                SEI SAFETY OFF
                    SORT CN Entry TO "Lookup entry.dbf"
SET SAFETY ON
USE "Lookup entry.dbf"
                                                                                                                                                          and the state of t
                              LOCATE FOR Look=Bobject
                              IF .. NOT. FOUND()
        CI FAR
LOOP
                            PULT
                             BROWSE
                    STORE Entry TO Searchval CLOSE DATABASES
                            ERASE . Lookup entry . dbf"
                      ENDIF"
                           IP Dobjecto'
SET EXACT OFF
SET SAFETY OFF
                            SORT ON descriptor TO *Loology descriptor.dbf*
                           SET SAFETY On
USE 'Lookup descriptor.dbf'
                           LOCATE FOR UPPER(TRIM(descriptor))=UPPER(TRIM(Dobject))
                           IF .NOT.FOUND()
```

```
1007
                                 RADIF
                                 BROWSE
                                STORE Entry TO Searchval CLOSE DATABASES
                                ERASE 'Lookup descriptor.dbf'
                                BET EXACT ON
                               ENDIP .
     CLEAR
                     ? 'Northern enalysis for entry '
       ?? Searchval
     ? 'Enter Y to proceed'
WAIT TO OK
CLEAR
IF UPPER (OK) <> 'Y'
sorreen 1 off
       screen 1 off
RETURN
RETURN
BNDIP

"COMPRESSION SUBROUTINE FOR Library.dbf
? 'Compressing the Libraries file now...'
USE 'SMARTGUY:FOXBASE+/Mac:Pox files:libraries.dbf'
SET SAFETY ON
USE 'Compressed libraries.dbf'
DELSTE FOR entered=0
PACK
                                                                            AND A STATE OF A SECOND STATE AND SECOND STATES AND AND SECOND SE
                             MARK1 = 1
                             5W2=0.
                            DO WHILE SW2=0 ROLL
                                 IF MARKI >= TOT
                               PACK
                           5W2=1
LCOP
                                  ENDIF
                            GO MARKI.
                       STORE library TO TESTA
                             STORE Library TO TESTS
                            IF TESTA = TESTB
DELETE
                            endif
                           MARK1 = MARK1+1 '
                            ENDDO ROLL
                       * Northern analysis
                           CLEAR
7 'Doing the northern now...
SET TALK ON
USE 'SmartGuy:FoxBASE+/Mac:Fox files:clones.dbf'
                           SET SAFETY OFF
COPY TO 'Hits.dbf' FOR entry-searchval
                             SET SAPETY ON
```

```
* MASTER ANALYSIS 3; VERSION 12-9-94
      Master menu for analysis output | 2005 falls of a saving to be on a to the saving saving and
    CLOSE DATABASES
    SET TALK OFF
    SET SAFETY OFF
    CLEAR
   SET DEVICE TO SCREEN
    SET DEFAULT TO "SmartGuy: FoxBASE+/Mac: fox files:Output programs:"
   USE "SmartGuy: FoxBASE+/Mac; fox, files: Clones, dbf" ... ) HERVE DE CAREER CONTROL OF THE CONTR
   STORE NUMBER TO INITIATE
    GO BOTTOM
    STORE NUMBER TO TERMINATE
  STORE 0 TO ENTIRE
STORE 0 TO CONDEN
STORE 0 TO ANAL.
   STORE O TO EMATCH
   STORE O TO HMATCH
   STORE 0 TO OMATCH
  STORE 0 TO IMATCH
STORE 0 TO XMATCH
   STORE O TO PRINTON
   STORE 0 TO PTF
  DO WHILE .T.
    Program.: Master analysis.fmt
   * Date....: 12/ 9/94
  * Version.: FoxBASE+/Mac, revision 1.10
  * Notes...: Format file Master analysis
SCREEN 1 TYPE 0 HEADING "Screen 1" AT 40.2 SIZE 286,492 PIXELS FONT "Geneva",9 COLOR 0.0.0. 6 PIXELS 39,255 TO 277,430 STYLE 28447 COLOR 0.0.-1,-25600,-1,-1 PIXELS 75,120 TO 178,241 STYLE 3871 COLOR 0.0.-1,-25600,-1,-1 PIXELS 27,98 SAY "Customized Output Menu" STYLE 65536 FONT "Geneva",274 COLOR 0.0.-1,-1,-1 PIXELS 45,54 GET condem STYLE 65536 FONT "Chicago",12 PICTURE "G*C Condensed format SIZE PIXELS 45,261 GET anal STYLE 65536 FONT "Chicago",12 PICTURE "G*C Condensed format SIZE PIXELS 54,261 GET anal STYLE 65536 FONT "Chicago",12 PICTURE "G*C Evact " SIZE 15.62 CO
PIXELS 117.126 GET EMATCH STYLE 65536 FONT "Chicago",12 FICTURE "8*C Exact " SIZE 15,62 CO PIXELS 135,126 GET HMATCH STYLE 65536 FONT "Chicago",12 FICTURE "8*C Homologous" SIZE 15,1 PIXELS 153,126 GET OMATCH STYLE 65536 FONT "Chicago",12 FICTURE "8*C Other spc" SIZE 15,84 PIXELS 90,152 SAY "Matches: " STYLE 65536 FONT "Geneva",268 COLOR 0,0,-1,-1,-1,-1
PIXELS 90,152 SAY "Matches: STYLE 65536 FONT "Geneva",268 COLOR 0,0,-1,-1,-1,-1

9 FIXELS 63,54 GET FRINTON STYLE 65536 FONT "Chicago",12 PICTURE "0*C Include clone listing"

6 PIXELS 171,126 GET Imatch STYLE 65536 FONT "Chicago",12 PICTURE "0*C Incyte" SIZE 15,65 CO

7 PIXELS 252,146 GET initiate STYLE 0 FONT "Geneva",12 SIZE 15,70 COLOR 0,0,-1,-1,-1,-1

9 PIXELS 270,146 GET terminate STYLE 0 FONT "Geneva",12 SIZE 15,70 COLOR 0,0,-1,-1,-1,-1

6 PIXELS 234,134 SAY "Include clones " STYLE 65536 FONT "Geneva",12 COLOR 0,0,-1,-1,-1,-1
@ PIXELS 270,125 SAY "->" STYLE 65536 FONT "Geneva",14 COLOR 0,0,-1,-1,-1,-1
@ PIXELS 198,126 GET PIF STYLE 65536 FONT "Chicago",12 PICTURE "@*C Print to file" SIZE 15,9
@ PIXELS 189,0 TO 257,120 STYLE 3871 COLOR 0,0,-1,-25600,-1,-1
@ PIXELS 209,8 SAY "Library selection" STYLE 65536 FONT "Geneva", 266 COLOR 0,0,-1,-1,-1,-1 @ PIXELS 227,18 GET ENTIRE STYLE 65536 FONT "Chicago", 12 PICTURE "@*RV All; Selected" SIZE 16
 * EOF: Master analysis.fmt
READ
      IF ANAL=9
      CLEAR
      CLOSE DATABASES
      ERASE TEMPMASTER.DBF
      USE "SmartGuy:FoxBASE+/Mac:fox files:clones.dbf"
      SET SAFETY ON
      SCREEN 1 OFF
      RETURN
      ENDIF
clear
? INITIATE
? TERMINATE
? CONDEN
? ANAL
```

```
? ematch
? Hmatch
? Cmatch
? IMATCH
SET TALK ON
           IF ENTIRES2
                                                                                                                                                                                          We market & March
     USE 'Unique libraries.dbf'
REPLACE ALL i WITH
BROWSE FIELDS i, libname, library, total, entered AT 0,0
     ENDIP
      USE "SmartGuy:FoxBASE+/Mac:fox files:clones.dbf"
      *COPY TO TEMPNUM FOR NUMBER>=INITIATE.AND.NAMER<=TERMINATE
    *USE TEMPNUM
    COPY STRUCTURE TO TEMPLIB
     USE TEMPLIB
           IF ENTIRE=1
   APPEND FROM 'SmartGuy: FoxBASE+/Mac:fox files:Clones.dbf'
ENDIF
                                                                                                                           The Section of the Se
            IF ENTIRE=2
    USE "Unique libraries.dbf"
           COPY TO SELECTED FOR UPPER(1)='Y' AND THE RESERVE OF THE SECOND CONTROL OF THE PROPERTY OF THE
            STORE RECCOUNT() TO STOPIT
         MARK=1
     DO WHILE .T.
                   IF MARK>STOPIT
                  CLEAR
              EXIT
                                                                                                                            ALTO ALL MODEL OF MATHEMATICAL ALL AND CALLS
           ENDIF
                 USE SELECTED
                  GO MARK
   STORE library TO THISONE
                   ? 'COPYING '
                  ?? THISONE
                  USE TEMPLIE
                  APPEND FROM "SmartGuy: FoxBASE+/Mac:fox files:Clones.dbf" FOR library=THISONE
                  STORE MARK+1 TO MARK
                 LOOP
                ENDDO
      ENDIP
  USE "SmartGuy:FoxBASE+/Mac:fox files:clones.dbf"
 COUNT TO STARTOT
 COPY STRUCTURE TO TEMPDESIG
 USE TEMPDESIG
  IF Ematch=0 .AND. Hmatch=0 .AND. Cmatch=0 .AND. IMATCH=0
       APPEND FROM TEMPLIB
 ENDIF
         IF Ematch=1
         APPEND FROM TEMPLIB FOR D='E'
        ENDIF
         IF Hmatchel
         APPEND FROM TEMPLIE FOR D='H'
        ENDIF
         IF Cmatchel
        APPEND FROM TEMPLIB FOR Da'O'
        ENDIP
       IF Imatch=1
        APPEND FROM TEMPLIB FOR D='I'.OR.D='X'.OR.D='N'
       ENDIF
        IF Xmatchel
       APPEND FROM TEMPLIB FOR D='X'
       ENDIF
COUNT TO ANALITOT
set talk off
DO CASE
```

```
CASE (PTP=0) [96] (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) (1.0) 
          SET DEVICE TO PRINT
          SET PRINT ON SAME SAME SAME SAME
          EJECT
         CASE PTF=1
         SET ALTERNATE TO "Total function sort.txt"
          *SET ALTERNATE TO "H and O function sort.txt"
         *SET ALTERNATE TO "Shear Stress HUVEC 2: Abundance sort.txt"
*SET ALTERNATE TO "Shear Stress HUVEC 2: Abundance con.txt"
         *SET ALTERNATE TO "Shear Stress HOVEC 2: Function sort.txt"
        *SET ALTERNATE TO "Shear Stress HUVEC 2:Distribution sort.txt"
*SET ALTERNATE TO "Shear stress HUVEC 1:Clone list.txt"
         *SET ALTERNATE TO "Shear Stress HUVEC 2:Location Bort.txt"
      SET ALTERNATE ON CALLERY.

ENDCASE CARROLLERY OF THE PROPERTY AND THE ALTERNATION OF THE PROPERTY OF THE PROPE
       IF PRINTON=1
       61,30 SAY "Database Subset Analysis" STYLE 65536 FONT "Geneva",274 COLOR 0,0,0,-1,-1,-1
       ENDIF WAY LEADING
       34 9 E
      ?:
                                                                                                                                                                               PROPERTY AND PARTY OF WILLIAMS FOR THE
     ? date() approved the larger of the province of the larger of the larger
     7 'Clone numbers '
      ?? STR(INITIATE, 6, 0)
     ?? 'through '
     ?? STR (TERMINATE, 6, 0)
      ? 'Libraries: !
     IF ENTIRE=1
                                                                                                          Bether the profite a profit of the control of the profit o
     ? 'All libraries'
    ENDIF
                                                                                 Commence of the second of
    IF ENTIRE=2
        MARK-1
           DO WHILE .T.
                        IF MARK-STOPIT
                        EXIT
                                                                                                                                   Section of the section of the
          ENDIF
USE SELECTED
                       ENDIF
       GO MARK
            ?? TRIM(libname)
                        STORE MARK+1 TO MARK
                        LOOP .
                                                                                   ENDDO
 ENDIF
  ? 'Designations: '
  IF Ematch=0 .AND. Hmatch=0 .AND. Omatch=0 .AND. IMATCH=0
 ?? 'All'
ENDIP
 IF Ematch=1
 ?? 'Exact,'
ENDIP
 IF Hmatch=1
 ?? 'Human, '
ENDIF .
IF Omatch=1
 ?? 'Other sp. '
EVDIF
IF Imatch=1
77 'INCYTE'
ENDIF
IF Xmatch=1
?? 'EST'
```

```
ENDIF
   ? 'Condensed format analysis'
                                              ENDIF
   IF ANALe1
  ? 'Sorted by NUMBER'
  ENDIF
   IF ANAL=2
  ? 'Sorted by ENTRY'
  IF ANALES
  ENDIF
                            great in the same and the
  IF ANAL=4
  ? 'Sorted by INTEREST'
 ENDIF
  IF ANAL=5
  7. !Arranged by LOCATION! He wis to have store participation of the second of the representation of the repres
 ENDIF
  IF ANAL=6
  ? 'Arranged by DISTRIBUTION'
 ENDIF
 IF ANAL=7
  ? 'Arranged by FUNCTION'
                                           •
 ENDIF
                                                                                                                                The second of the contract of the second of the second of
                                                                                                    No. of Belleville
 ? 'Total clones represented: '
 ?? STR(STARTOT, 6, 0)
 ? 'Total clones analyzed: '
 ?? STR(ANALTOT, 6, 0)
 7 'l' = library d = designation f = distribution z = location r = function
 ********
 USE TEMPDESIG
 SCREEN 1 TYPE 0 HEADING "Screen 1" AT 40,2 SIZE 286,492 PIXELS FONT "Geneva",7 COLOR 0,0,0,
DO CASE
CASE ANAL=1
 * sort/number
SET HEADING ON
                                                                                              and the second of the second o
IF CONDEN=1
SORT TO TEMP1 ON ENTRY, NUMBER
DO "COMPRESSION number.PRG"
ELSE .
SORT TO TEMP1 ON NUMBER
USE TEMP1
list off fields number, L, D, F, Z, R, C, EMTRY, S, DESCRIPTOR
 *list off fields number, L, D, F, Z, R, C, ENTRY, S, DESCRIPTOR, LENGTH, RFEND, INIT, I
CLOSE DATABASES
ERASE TEMP1.DBF
ENDIF
CASE ANAL=2
 * sort/DESCRIPTOR
SET HEADING ON
*SORT TO TEMP1 ON DESCRIPTOR, ENTRY, NUMBER/S for D='E'.OR.D='H'.OR.D='O'.OR.D='X'.OR.D='I'
*SORT TO TEMP1 ON ENTRY, DESCRIPTOR, NUMBER/S for D='E'.OR.D='H'.OR.D='O'.OR.D='X'.OR.D='I'
SORT TO TEMP1 ON ENTRY, START/S for D='E'.OR.D='H'.OR.D='O'.OR.D='X'.OR.D='I'
IF CONDEN-1
DO "COMPRESSION entry PRG"
ELSE
  list off fields number, L, D, F, Z, R, C, ENTRY, S, DESCRIPTOR, LEXGTH, RFEND, INIT, I
CLOSE DATABASES
ERASE TEMP1.DAF
ENDIF
```

```
CASE ANAL=3 ANAL=3 ANAL=3
   * sort by abundance
   SET HEADING ON ...
   SORT TO TEMP1 ON ENTRY, NUMBER for D='E4 OR.D='H', OR.D='O'.OR.D='X', OR.D='I'
  DO "COMPRESSION abundance.PRG"
  CASE ANAL-4
   f sort/interest
   SET HEADING ON
  IF CONDEN=1
  SORT TO TEMP1 ON ENTRY, NUMBER FOR I>0
  DO "COMPRESSION interest.PRG"
  elşe
  SORT ON I/D, ENTRY TO TEMP1 FOR I>1
  USB TEMP1
  list off fields number, L, D, F, Z, R, C, ENTRY, S, DESCRIPTOR, LENGTH, RFEND, INIT, I
  CLOSE DATABASES
  ERASE TEMP1.DBF
  ENDIF
  CASE ANAL=5...
                                                                            A THE RESERVE OF LICENSESSES THE PROPERTY OF THE PARTY OF
  * arrange/location
  SET HEADING ON
                                                                                                                       「ジン」
(4:8 ) - 伊藤女子
                                                                                                                                                      . . .
  STORE 4 TO AMPLIFIER
                                                                                                                                                            Procession Charles
                                                                                                       1. 1 B.
  ? 'Nuclear:'
  SORT ON ENTRY, NUMBER FIELDS REEND, NUMBER, L.D. F.Z.R, C. ENTRY, S. DESCRIPTOR, LENGTH, INIT, I, COMMEN
 IF CONDEN=1
 DO. "Compression location.prg"
 ELSE
 DO "Normal subroutine 1"
 ENDIP
 ? 'Cytoplasmic:'
 SORT ON ENTRY, NUMBER FIELDS REEND, NUMBER, L. D. F. Z. R. C. ENTRY, S. DESCRIPTOR, LENGTH, INIT, I, COMMEN
 IF CONDEN=1
 DO "Compression location.prg"
 BLSE
 DO "Normal subroutine 1"
PATITR
 ? 'Cytoskeleton:'
 SORT ON ENTRY, NUMBER FIELDS RFEND, NUMBER, L.D.F.Z.R.C. ENTRY, S. DESCRIPTOR, LENGTH, INIT, I. COMMEN
 IF CONDENS:
DO 'Compression location.prg'
ELSE
DO "Normal subroutine 1"
ENDTE
? 'Cell surface: '
SORT ON ENTRY, NUMBER FIELDS REEND, NUMBER, L, D, F, Z, R, C, ENTRY, S, DESCRIPTOR, LENGTH, INIT, I, COMMEN
IF CONDEN=1
DO "Compression location.prg"
ELSE
DO "Normal subroutine 1"
ENDIF
? 'Intracellular membrane:'
SORT ON ENTRY, NUMBER FIELDS RFEND, NUMBER, L.D. F. Z.R.C. ENTRY, S. DESCRIPTOR, LENGTH, INIT, I, COMMEN
IF CONDEN=1
DO "Compression location.prg"
ELSE
DO "Normal subroutine 1"
PNDIP
? 'Mitochondrial:'
SORT ON ENTRY, NUMBER FIELDS RFEND, NUMBER, L, D, F, Z, R, C, ENTRY, S, DESCRIPTOR, LENGTH, INIT, I, COMMEN
IF CONDEN-1
DO "Compression location.prg"
ELSE.
DO "Normal subroutine 1"
ENDIP
```

```
? Secreted? Stew Time Fit Wellson Buckeys
 SORT ON ENTRY, NUMBER FIELDS RFEND, NUMBER, L.D.F.Z.R.C. ENTRY, S. DESCRIPTOR, LENGTH, INIT, I, COMMEN
IF CONDEN=1
DO "Compression location.prg"
 ELSE
 DO "Normal subroutine 1"
ENDIF
 ? !Other!
 SORT ON ENTRY, NUMBER FIELDS RFEND, NUMBER, L. D. F. Z. R. C. ENTRY, S. DESCRIPTOR, LENGTH, INIT, I, COMMEN
IF CONDEN=1
DO "Compression location.prg"
 ELSE
 DO "Normal subroutine 1"
ENDIF:
 ? 'Unknown:'
SORT ON ENTRY, NUMBER FIELDS RFEND, NUMBER, L.D. F. Z.R.C. ENTRY, S. DESCRIPTOR, LENGTH, INIT, I, COMMEN
 IF CONDEN-1
DO "Compression location.prg"
ELSE
DO "Normal subroutine 1"
ENDIF
IF CONDEN-1
SET DEVICE. TO PRINTER
SET PRINTER ON
EJECT
DO "Output heading.prg"
USE Analysis location.dbf
DO "Create bargraph.prg"
SET HEADING OFF
          FUNCTIONAL CLASS
                                                  TOTAL
                                                          UNIQUE NEW $ TOTAL
LIST OFF FIELDS Z, NAME, CLONES, GENES, NEW, FERCENT, GRAPH
CLOSE DATABASES
ERASE TEMP2. DBF
SET HEADING ON
*USE *SmartGuy:FoxBASE+/Mac:fox files:TEMFMASTER.dbf*
ENDIF .
CASE ANAL=6
* arrange/distribution
SET HEADING ON
STORE 3 TO AMPLIFIER
? 'Cell/tissue specific distribution:'
SORT ON ENTRY, NUMBER FIELDS RFEND, NUMBER, L.D.F.Z.R.C. ENTRY, S. DESCRIPTOR, LENGTH, INIT, I, COMMEN
IF CONDENSI
DO "Compression distrib.prg"
ELSE
DO "Normal subroutine 1"
ENDIF
7 'Non-specific distribution:'
SORT ON ENTRY, NUMBER FIELDS RFEND, NUMBER, L.D.F.Z.R.C. ENTRY, S. DESCRIPTOR, LENGTH, INIT, I. COMMEN
IF CONDEN=1
DO "Compression distrib.prg"
ELSE
DO "Normal subroutine 1"
ENDIF
? 'Unknown distribution:'
SORT ON ENTRY, NUMBER FIELDS RFEND, NUMBER, L.D. F.Z.R.C. ENTRY, S. DESCRIPTOR, LENGTH, INIT, I, COMMEN
IF CONDEN=1
DO "Compression distrib.prg"
ELSE
DO "Normal subroutine 1"
ENDIF
IF CONDEN=1
SET DEVICE TO PRINTER
SET PRINTER ON
```

```
EJECT - Grant Reading. Pro Control of Manager DO Coutput heading.prg*
USE "Analysis distribution.dbf*
 DO: Create bargraph.prg
  SET HEADING OFF
            FUNCTIONAL CLASS
                                                        TOTAL
                                                                 UNIQUE & TOTAL
 LIST OFF PIELDS P, NAME, CLONES, GENES, PERCENT, GRAPH
 CLOSE DATABASES
 ERASE TEMP2 DBF
 SET HEADING ON
  *USE *SmartGuy: FoxBASE+/Mac:fox files: TEMPMASTER.dbf*
 ENDIF
 CASE ANAL=7
 * arrange/function
 SET HEADING ON
 STORE 10 TO AMPLIFIER
                                     BINDING PROTEINS
 ? 'Surface molecules and receptors:'
SORT ON ENTRY, NUMBER FIELDS RFEND, NUMBER, L, D, F, Z, R, C, ENTRY, S, DESCRIPTOR, LENGTH, INIT, I, COMMEN
 DO *Compression function.prg*
 ELSE
 DO "Normal subroutine 1"
 ENDIF
 ? 'Calcium-binding proteins:'
 SORT ON ENTRY, NUMBER FIELDS RFEND, NUMBER, L, D, F, Z, R, C, ENTRY, S, DESCRIPTOR, LEXGIH, INIT, I, COMMEN
 IF CONDEN=1
 DO *Compression function.prg*
 ELSE
 DO 'Normal subroutine 1"
ENDIF
 ? 'Ligands and effectors:'
SORT ON ENTRY NUMBER FIELDS REEND, NUMBER, L, D, F, Z, R, C, ENTRY, S, DESCRIPTOR, LEXSTH, INIT, I, COMMEN
IF CONDEN=1
DO 'Compression function.prg'
ELSE
DO 'Normal subroutine 1'
ENDIF
? 'Other binding proteins:'
SORT ON ENTRY, NUMBER FIELDS RFEND, NUMBER, L.D. F, Z, R, C, ENTRY, S, DESCRIPTOR, LENGTH, INIT, I, COMMEN
IP CONDEN=1
DO "Compression function.prg"
ELSE
DO "Normal subroutine 1"
ENDIF
*EJECT
?-1,_
                                    ONCOGENES!
? 'General oncogenes:'
SORT ON ENTRY, NUMBER FIELDS RFEND, NUMBER, L, D, F, Z, R, C, ENTRY, S, DESCRIPTOR, LENGTH, INIT, I, COMMEN
IF CONDEN=1
DO "Compression function.prg"
ELSE
DO 'Normal subroutine 1"
ENDIF
7 'GTP-binding proteins:'
SORT ON ENTRY, NUMBER FIELDS RPEND, NUMBER, L, D, F, Z, R, C, ENTRY, S, DESCRIPTOR, LENGTH, INIT, I, COMMEN
IF CONDEN=1
DO 'Compression function.prg'
ELSE
DO 'Normal subroutine 1"
ENDIP
? 'Viral elements:'
```

```
SORT ON ENTRY, NUMBER FIELDS RFEND, NUMBER, L, D, P, Z, R, C, ENTRY, S, DESCRIPTOR, LENGTH, INIT, I, COMMEN IF CONDEN-1
 ELSE
 DO "Normal subroutine 1"
 ENDIF
 ? 'Kinases and Phosphatases:'
 SORT ON ENTRY, NUMBER FIELDS REEND, NUMBER, L.D. F. Z. R. C. ENTRY, S. DESCRIPTOR, LENGTH, INIT, I, COMMEN
 IF CONDEN=1
 Do "Compression function.prg"
 DO "Normal subroutine 1"
 ENDIF
 ? 'Tumor-related antigens:'
 SORT ON ENTRY, NUMBER FIELDS RFEND, NUMBER, L.D. F.Z. R.C. ENTRY, S, DESCRIPTOR, LENGTH, INIT, I, COMMEN
 IF CONDEN=1
 DO "Compression function.prg"
 ELSE
DO "Normal subroutine 1"
ENDIF
 *EJECT
 7 '
                              PROTEIN SYNTHETIC MACHINERY PROTEINS!
? 'Transcription and Nucleic Acid-binding proteins:'
SORT ON ENTRY, NUMBER FIELDS RFEND, NUMBER, L, D, F, Z, R, C, ENTRY, S, DESCRIPTOR, LENGTH, INIT, I, COMMEN
IF CONDEN-1
DO "Compression function.prg"
ELSE
DO 'Normal subroutine 1'
ENDIF
? 'Translation:'
SORT ON ENTRY, NUMBER FIELDS RFEND, NUMBER, L, D, F, Z, R, C, ENTRY, S, DESCRIPTOR, LENGTH, INIT, I, COMMEN
IF CONDEN=1
DO "Compression function.prg"
ELSE
DO 'Normal subroutine 1'
ENDIF
? 'Ribosomal proteins:'
SORT ON ENTRY, NUMBER FIELDS REEND, NUMBER, L, D, F, Z, R, C, ENTRY, S, DESCRIPTOR, LENGTH, INIT, I, COMMEN
IF CONDENS1
DO "Compression function.prg"
ELSE
DO 'Normal subroutine 1'
ENDIF
? 'Protein processing:'
sort on entry, number fields rpend, number, L, D, F, Z, R, C, Entry, S, Descriptor, Length, Init, I, commen
IF CONDEN=1
DO "Compression function.prg".
ELSE
DO 'Normal subroutine l'
ENDIF
*EJECT
                                ENZYMES'
? 'Ferroproteins:'
SORT ON ENTRY, NUMBER FIELDS REEND, NUMBER, L. D. F. Z. R. C. ENTRY, S. DESCRIPTOR, LENGTH, INIT, I, COMMEN
IF CONDEN=1
DO "Compression function.prg"
ELSE
DO 'Normal subroutine 1'
ENDIF
? 'Proteases and inhibitors:'
SORT ON ENTRY, NUMBER FIELDS RFEND, NUMBER, L.D.F.Z.R.C. ENTRY, S.DESCRIPTOR, LENGTH, INIT, I, COMMEN
IF CONDEN=1
DO "Compression function.prg"
ELSE
```

```
ENDIP
 "? "Oxidative phosphorylation:
  SORT ON ENTRY, NUMBER FIELDS RFEND, NUMBER, L, D, F, Z, R, C, ENTRY, S, DESCRIPTOR, LENGTH, INIT, I, COMMEN
  IF CONDEN=1
  DO "Compression function.prg"
  ELSE -
  DO "Normal subroutine 1"
  ENDIP
  7 Sugar metabolism:
  SORT ON ENTRY, NUMBER FIELDS REEND, NUMBER, L, D, F, Z, R, C, ENTRY, S, DESCRIPTOR, LENGTH, INIT, I, COMMEN
  IF CONDEN-1
  DO "Compression function.prg"
 ELSZ
 DO Normal subroutine 1"
 ENDIF
  ? 'Amino acid metabolism:'
 SORT ON ENTRY, NUMBER FIELDS RFEND, NUMBER, L. D. F. Z. R. C. ENTRY, S. DESCRIPTOR, LENGTH, INIT, I, COMMEN
 IF CONDEN=1
 DO "Compression function.prg"
 EZSE
 DO "Normal subroutine 1"
 ENDIF
 ? 'Nucleic acid metabolism:'
 SORT ON ENTRY, NUMBER FIELDS REEND, NUMBER, L.D. F. Z. R, C, ENTRY, S, DESCRIPTOR, LENGTH, INIT, I, COMMEN
 IF CONDEN=1
 DO *Compression function.prg*
 ELSE
 DO 'Normal subroutine 1'
 ENDIP
? 'Lipid metabolism:'
 SORT ON ENTRY, NUMBER FIELDS RFEND, NUMBER, L, D, F, Z, R, C, ENTRY, S, DESCRIPTOR, LENGTH, INIT, I, COMMEN
 IF CONDEN=1
DO "Compression function.prg"
 ELSE:
 DO 'Normal subroutine 1'
 ENDIP
 ? 'Other enzymes:'
SORT ON ENTRY, NUMBER FIELDS REEND, NUMBER, L. D. F. Z. R. C. ENTRY, S. DESCRIPTOR, LENGTH, INIT, I, COMMEN
IF CONDEN=1
DO "Compression function.prg"
 ELSE
DO "Normal subroutine 1"
ENDIF
 *EJECT
?
                                 MISCELLANEOUS CATEGORIES'
7 'Stress'response:'
SORT ON ENTRY, NUMBER FIELDS REEND, MIMBER, L. D. F. Z. R. C. ENTRY, S. DESCRIPTOR, LENGTH, INIT, I. COMMEN
IF CONDEN=1
DO "Compression function.prg"
ELSE
DO 'Normal subroutine 1"
ENDIF
7 'Structural:'
SORT ON ENTRY, NUMBER FIELDS REEND, NUMBER, L, D, F, Z, R, C, ENTRY, S, DESCRIPTOR, LENGTH, INIT, I, COMMEN
IF COMDEN=1
DO "Compression function.prg"
ELSE
DO "Normal subroutine 1"
ENDIF
? 'Other clones:'
SORT ON ENTRY NUMBER FIELDS REEND, NUMBER, L, D, P, Z, R, C, ENTRY, S; DESCRIPTOR, LENGTH, INIT, I, COMMEN
IF CONDEN=1
DO *Compression function.prg*
ELSE
```

```
-DO: Normal subroutine 1 18 3 17 25 8 1 - 4 Collection
PENDIF
(?//Clones of unknown function:
 SORT ON ENTRY, NUMBER FIELDS REEND, NUMBER, L.D.F.Z.R.C.ENTRY, S, DESCRIPTOR, LENGTH, INIT, I, COMMEN
IF CONDEN=1
DO Compression function.prg*
ELSE
DO "Normal subroutine 1"
 ENDIF.
IF CONDENSATION AND DO NOT BE THE
EJECT
 *SET DEVICE TO PRINTER
*SET PRINT ON
DO "Output heading.prg"
USE "Analysis function.dbf"
DO "Create bargraph.prg"
SET HEADING OFF
***
SCREEN 1 TYPE 0 HEADING "Screen 1" AT 40,2 SIZE 286,492 PIXELS FONT "Geneva",12 COLOR 0,0,0
    4.
2 .
                                                                 TOTAL
                                                                          TOTAL
                                                                                   NEW
                                                                                             DIST
         FUNCTIONAL CLASS
                                                   CLONES
                                                            GENES GENES
                                                                             FUNCTIONAL CLASS'
7 1/
*LIST OFF FIELDS P, NAME, CLONES, GENES, NEW, PERCENT, GRAPH, COMPANY
LIST OFF FIELDS P, NAME, CLONES, GENES, NEW, PERCENT, GRAPH
CLOSE DATABASES
ERASE TEMP2.DBF
SET HEADING ON
*USE *SmartGuy:FoxBASE+/Mac:fox files:TEMPMASTER.dbf*
ENDIF
CASE ANAL 8
DO "Subgroup summary 3.prg"
ENDCASE
DO "Test print.prg"
SET FRINT OFF
SET DEVICE TO SCREEN
CLOSE DATABASES
*ERASE TEMPLIB.DBF
*ERASE TEMPNUM.DBF
*ERASE TEMPDESIG.DEF
*ERASE SELECTED.DBF
CLEAR
LCOP
ENDDO
```

```
* COMPRESSION SUBROUTINE FOR ANALYSIS PROGRAMS
 USE TEMP1
 COUNT TO TOT
 REPLACE ALL REEND WITH 1
 MARK1 = 1
 5W2=0
 DO WHILE SW2=0 ROLL
   IF MARK1 >= TOT
   PACK
   COUNT TO UNIQUE
COUNT TO NEWGENES FOR D='H'.OR.D='O'
   .$W2=1
   LOOP
   ENDIF
GO MARKI
DUP = 1
STORE ENTRY TO TESTA
5W = 0
DO WHILE SW=0 TEST
SKIP
STORE ENTRY TO TESTE
  IF TESTA = TESTB
  DELETE
  DUP = DUP+1
  LOOP
  ENDIF
GO MARKI
REPLACE REEND WITH DUP
MARK1 = MARK1+DUP
SN=1
LOOP
ENDDO TEST
LOOP
ENDDO ROLL
GO TOP
STORE Z TO LOC
USE 'Analysis location.dbf' LOCATE FOR Z-LOC
REPLACE CLONES WITH TOT
REPLACE GENES WITH UNIQUE
REPLACE NEW WITH NEWGENES
USE TEMP1
SORT ON RFEND/D TO TEMP2
USE TEMP2
?? STR(UNIQUE, 5, 0)
?? genes, for a total of '?? STR(TOT,5,0)
??'.clones'
                            V Coincidence'
list off fields number, RFEND, L. D. F. Z. R. C. ENTRY, S. DESCRIPTOR, LENGTH, INIT, I
*SET PRINT OFF
CLOSE DATABASES
ERASE TEMP1.DBF
ERASE TEMP2.DBF
USE TEMPDESIG
```

```
* COMPRESSION SUBROUTINE FOR ANALYSIS PROGRAMS
 USE TEMP1
 COUNT TO TOT
 REPLACE ALL RFEND WITH 1
 MARK1 = 1
 5W2=0
 DO WHILE SW2=0 ROLL
   IF MARK1 >= TOT
   PACK
   COUNT TO UNIQUE
   5W2=1
   LOOP
   ENDIF
 GO MARKI
 DUP = 1
 STORE ENTRY TO TESTA
 SW = 0
 DO WHILE SW=0 TEST
 SKIP
 STORE ENTRY TO TESTE
   IF TESTA = TESTB
   DELETE
   DUP = DUP+1
   LOOP
  ENDIF
GO MARKI
REPLACE RFEND WITH DUP
MARK1 = MARK1+DUP
SW-1
LOOP
ENDDO TEST
LOOP
ENDDO ROLL
*BROWSE
*SET PRINTER ON
SORT ON DATE TO TEMP2
USE TEMP2
?? STR(UNIQUE, 4,0)
?? ' genes, for a total of'
?? STR(TOT,4,0)
?? " clones"
                           V Coincidence'
COUNT TO P4 FOR I=4
IF P4>0
? STR(P4,3,0)
?? ' genes with priority = 4 (Secondary analysis:)'
list off fields number,RFEND,L,D,F,Z,R,C,ENGRY,S,DESCRIPTOR,LENGTH,INIT for I=4
ENDIF
COUNT TO P3 FOR I=3
IF P3>0
? STR(P3,3,0)
?? ' genes with priority = 3 (Full insert sequence:)'
list off fields number.RFEND, L.D, F, Z, R, C, EMTRY, S, DESCRIPTOR, LEXCTH, INIT for I=3
ENDIP
COUNT TO P2 FOR I=2.
IF P2>0
? STR(P2,3,0)
?? ' genes with priority = 2 (Primary analysis complete:)'
list off fields number, RFEND, L, D, F, Z, R, C, EMTRY, S, DESCRIPTOR, LENGTH, INIT for I=2
ENDIF
COUNT TO PL FOR Isl
IF P1>0
```

```
7.STR(P1,3,0)
?? genes with priority = 1 (Primary snallysis needed:)'
List off fields number.RFEND, L, D, P, Z, R, C, ENTRY, S, DESCRIPTOR, LENGTH, INIT for I=1
              ENDIF
           *SET PRINT OFF CONTROL OF THE CONTRO
              USE 'SmartGuy:FoxBASE+/Mac:fox files:clones.dbf'
                 And Silver of the Silver of th
                                                    **
                 5017
        . .
                    41.4
        12 1
A CONTRACTOR OF THE STATE OF TH
```

```
* COMPRESSION SUBROUTINE FOR ANALYSIS PROGRAMS...
     USE: TEMP1:
     COUNT TO TOTAL AND THE STATE OF THE COUNTY O
     MARK1 = 1
     5W2=0
     DO WHILE SW2=0 ROLL
            IF MARKL >= TOT
           PACK
           COUNT TO UNIQUE ;
           5W2=1
           LOOP
          ENDIP IN THE PROPERTY.
    GO MARKI
    DUP = 1
    STORE ENTRY TO TESTA
    SW = 0 ...
    DO WHILE 5W=0 TEST
   SKIP
    STORE ENTRY TO TESTS
          IF TESTA = TESTB
          DELETE
         DUP = DUP+1
LOOP
ENDIF
   GO MARKI
  REPLACE RFEND WITH DUP
  MARK1 = MARK1+DUP
  SW=1
 LOOP
ENDDO TEST
  LOOP
  ENDDO ROLL
  *BROWSE
  *SET PRINTER ON
  SORT ON NUMBER TO TEMP2
 USZ TEMP2
?? ' genes, for a total of '
?? STR(TOT,5,0)
?? ' clones'
                                                                          V Coincidence'
list off fields number.RFEND,L,D,F,Z,R,C,ENTRY,S,DESCRIPTOR,LENGTH,INIT,I
*SET PRINT OFF
CLOSE DATABASES
ERASE TEMP1.DBF
ERASE TEMP2.DBF
USE 'SmartGuy:FoxBASE+/Mac:fox files:clones.dbf'
```

```
* COMPRESSION SUBROUTINE FOR ANALYSIS PROGRAMS
USE TEMP1
 COUNT TO TOT
 REPLACE ALL REEND WITH 1
MARK1 = 1
5W2=0
COUNT TO UNIQUE
COUNT TO NEWGENES FOR D='H'.OR.D='O'
  6W2=1
  LOOP
  ENDIF
GO MARKI
DUP - 1
STORE ENTRY TO TESTA
SW = 0
DO WHILE SW=0 TEST
STORE ENTRY TO TESTE
  IF TESTA = TESTB
DELETE
  DUP = DUP+1
  LOOP
  ENDIF
GO MARKI
REPLACE RFEND WITH DUP
MARK1 = MARK1+DUP
SW=1
LOOP
ENDDO TEST
LOOP
ENDDO ROLL
GO TOP
STORE R TO FUNC
USE "Analysis function.dbf"
LOCATE FOR P=FUNC
REFLACE CLONES WITH TOT
REPLACE GENES WITH UNIQUE
REPLACE NEW WITH NEWGENES.
USE TEMP1
SORT ON RFEND/D TO TEMP2
USE TEMP2
SET HEADING ON
?? STR (UNIQUE, 5, 0)
?? ' genes, for a total of '
?? STR(TOT,5,0)
?? 'clones'
***
? '
                        V Coincidence'
list off fields number, RFEND, L, D, F, Z, R, C, EMTRY, S, DESCRIPTOR, LENGTH, INIT, I
*SCREEN 1 TYPE 0 HEADING "Screen 1" AT 40,2 SIZE 286,492 PIXELS FONT "Geneva",12 COLOR 0,0,
*list off fields RFEND, S, DESCRIPTOR
*SET PRINT OFF
CLOSE DATABASES
ERASE TEMP1.DBF
ERASE TEMP2.DBF
USE TEMPDESIG
```

```
* COMPRESSION SUBROUTINE FOR ANALYSIS PROGRAMS
    USB TEMP1
     COUNT TO TOT
    REPLACE ALL REEND WITH 1
   MARK1 = 1 ... ... 100 february recorded uniform the SW2=0. The least the least the DO WHILE SW2=0 ROLL
        IF MARK1 >= TOT
      PACK
COUNT TO UNIQUE
   SW2=1 Live to the server for LOOP 1975 to the server for L
DUF = 1
STORE ENTRY TO TESTA
SW = 0
DO WHILE SW=0 TEST
SKIP
STORE ENTRY TO TESTB
IF TESTA = TESTB
DELETE
DUP = DUP+1
LOOP
ENDIF
GO MARK1
REPLACE RFEND WITH DUP
MARK1 = MARK1+DUP
SW=1
  LOOP
  ENDDO TEST
  LOOP
  ENIDDO ROLL
  GO TOP .
  STORE F TO DIST
  USE "Analysis distribution.dbf"
 LOCATE FOR P-DIST
 REPLACE CLONES WITH TOT
 REPLACE GENES WITH UNIQUE
  USE TEMP1
 sort on rfend/d to TEMP2
 USE TEMP2
 ?? STR(UNIQUE,5,0)
?? 'genes, for a total of '
?? STR(TOT,5,0)
  ?? 'clones'
                                                                                                  V Coincidence'
 list off fields number, RFEND, L, D, F, Z, R, C, ENTRY, S, DESCRIPTOR, LENGTH, INIT, I
 *SET PRINT OFF
 CLOSE DATABASES
 ERASE TEMP1.DEF
ERASE TEMP2.DBF
 USE TEMPDESTO
```

```
COMPRESSION SUBROUTING FOR ANALYSIS PROGRAMS - TYPE TO BE ALL TO BE COMPRESSION SUBROUTING FOR ANALYSIS PROGRAMS - TYPE TO BE ALL TO BE COMPRESSION OF THE COMPRESSIO
       USE TEMP1
        COUNT TO TOT
       REPLACE ALL REEND WITH 1
       MARK1 = 1
       SW2=0
       DO WHILE SW2=0 ROLL
       / IF MARK1 >= TOT
              PACK
       COUNT TO UNIQUE
          SW2=1
       To LOOP (1994) And the control of the
            ENDIP
     GO MARKI
     DUP = 1
     STORE ENTRY TO TESTA
     0. = Wa
     DO WHILE SW=0 TEST
     SKIP.
     STORE ENTRY TO TESTE
           IF TESTA = TESTB
DELETE
            DUP = DUP+1
           LOOP
          ENDIF
   GO MARK1
  REPLACE REEND WITH DUP
MARK1 = MARK1+DUP
   SW=1
  LOOP .
  ENDDO TEST
  LOOP
  ENDDO ROLL
 GO TOP
 USB TEMP1
 7? STR(UNIQUE, 5,0)
77 ' genes, for a total of '
2? STR(TOT, 5,0)
 ??' clones'
                                                                                                     V Coincidence'
 list off fields number, RFEND, L, D, F, Z, R, C, ENTRY, S, DESCRIPTOR, LENGTH, INIT, I
 *SET PRINT OFF
CLOSE DATABASES
ERASE TEMP1.DBF
USE TEMPDESIG
```

```
COMPRESSION SUBROUTINE FOR ANALYSIS PROGRAMS
  USE SmartGuy: FoxBASE+/Mac: fox files: Clones.dbf*
COPY TO TEMP1 FOR
USE TEMP1
COUNT TO IDGENE FOR D='E'.OR.D='O'.OR.D='H'.OR.D='N'.OR.D='R'.OR.D='A'
  DELETE FOR D='N'.OR.D='D'.OR.D='A'.OR.D='U'.OR.D='S'.OR.D='M'.OR.D='R'.OR.D='V'
  PACK
  COUNT TO TOT
  REPLACE ALL REEND WITH 1
  MARK1 = 1
  57/2=0
  DO WHILE SW2=0 ROLL
    IF MARK1 >= TOT
    PACK
    COUNT TO UNIQUE
    SW2=1
    LOOP
    ENDIF
  GO MARKI
  DUP = 1
  STORE ENTRY TO TESTA
  5W = 0
 DO WHILE SW=0 TEST
 SKIP
 STORE ENTRY TO TESTE
    IF TESTA = TESTB
   DELETE.
   DUP = DUP+1
   LOOP
   ENDIF
 GO MARKI
 REPLACE RFEND WITH DUP
 MARK1 = MARK1+DUP
 5W=1
 LOOP
 ENDDO TEST
 LOOP
 ENDDO ROLL
 *BROWSE
 *SET PRINTER ON
 SORT ON REEND/D, NUMBER TO TEMP2
 USE TEMP2
 REPLACE ALL START WITH RFEND/IDGENE*10000
 ?? STR(UNIQUE,5,0)
 7? ' genes, for a total of '
7? STR(TOT,5,0)
?? ' clones'
? ' Coincidence V
                           V Clones/10000'
 set heading off
SCREEN 1 TYPE 0 HEADING "Screen 1" AT 40,2 SIZE 286,492 PIXELS FONT "Geneva",7 COLOR 0,0,0,
 list fields number, RFEND, START, L, D, F, Z, R, C, ENTRY, S, DESCRIPTOR, INIT, I
 *SET PRINT OFF
CLOSE DATABASES
ERASE TEMP1.DBF
ERASE TEMP2.DBF
USE *SmartGuy:FoxBASE+/Mac:fox files:clones.dbf*
```

```
COMPRESSION SUBROUTINE FOR ANALYSIS PROGRAMS
CUSE TEMP1 -
OCCUMT TO IDGENE FOR D='E'.OR.D='O'.OR.D='H'.OR.D='N'.OR.D='R'.OR.D='A'
DELETE FOR Da'N'.OR.D='D'.OR.D='A'.OR.D='U'.OR.D='S'.OR.D='M'.OR.D='R'.OR.D='V'
COUNT TO TOT
∴ SW2=0 \....
F DO WHILE SW2=0 ROLL
IF MARKI >= TOT
   PACK
   COUNT TO UNIQUE
   SW2=1
   LOOP
   ENDIF
 GO MARKI
 DUP = 1
 STORE ENTRY TO TESTA
 SW = 0
 DO WHILE SW=0 TEST
 SKIP
 STORE ENTRY TO TESTE
  IF TESTA = TESTE
DELETE
  DUP = DUP+1
  LOOP
  ENDIF
GO MARKI
REPLACE REEND WITH DUP
MARK1 = MARK1+DUP
5W=1
LOOP
ENDDO TEST
LOOP
ENDDO ROLL
*BROWSE
*SET PRINTER ON
SORT ON RFEND/D, NUMBER TO TEMP2
USE TEMP2
REPLACE ALL START WITH RFEND/IDGENE*10000
?? STR (UNIQUE, 5, 0)
?? ' genes, for a total of '
?? STR(TOT, 5, 0)
?? ' clones'
? ' Coincidence V
                          V Clones/10000'
set heading off
SCREEN 1 TYPE 0 HEADING "Screen 1" AT 40,2 SIZE 286,492 PIXELS FONT "Geneva",7 COLOR 0,0,0, list fields number, RFEND, START, L, D, F, Z; R, C, ENTRY, S, DESCRIPTOR; INIT, I
*SET PRINT OFF
CLOSE DATABASES
ERASE TEMP1.DBF
ERASE TEMP2.DBF
USE "SmartGuy: FoxBASE+/Mac: fox files:clones.dbf"
```

```
USE TEMP1
COUNT TO TOT
?? 'Total of'
?? STR(TOT,4,0)
?? 'clones'
?'
*list off fields number, L, D, F, Z, R, C, ENTRY, DESCRIPTOR, LENGTH, RFEND, INIT, I
list off fields number, L, D, F, Z, R, C, ENTRY, DESCRIPTOR
CLOSE DATABASES
ERASE TEMP1, DBF
USE TEMPDESIG
```

ENDDO

```
*Lifescan menu, version 8-7-94
  SET TALK OFP
  set device to screen
  CLEAR
  USE "SmartGuy: FoxBASE+/Mac:fox files:clones.dbf"
  STORE LUFDATE() TO Update
  GO BOTTOM
                                                                            All Deglet.
  STORE RECNO() TO cloneno
  STORE 6 TO Chooser
DO WHILE .T.
  * Program : Lifeseq menu.fmt
  * Date.... 1/11/95
 Version :: FoxEASE+/Mac, revision 1.10
Notes . . . :: Format file Lifeseg menu
SCREEN 1 TYPE 0 HEADING "Screen 1" AT 40,2 SIZE 286,492 PIXELS FONT "Geneva",268 COLOR 0,0,6 PIXELS 18,126 TO 77,365 STYLE 28479 COLOR 0,0,-1,-25600,-1,-16223,-16721,-15725 PIXELS 110,29 TO 188,217 STYLE 3871 COLOR 0,0,-1,-25600,-1,-1

PIXELS 45,161 SAY "LIFESEQ" STYLE 65536 FONT "Geneva",536 COLOR 0,0,-1,-1,7135,5884 PIXELS 36,269 SAY "TM" STYLE 65536 FONT "Geneva",12 COLOR 0,0,-1,-1,7135,5884 PIXELS 63,143 SAY "Molecular Biology Desktop" STYLE 65536 FONT "Helvetica",18 COLOR 0,0,0,0 PIXELS 90,252 TO 251,467 STYLE 28447 COLOR 0,0,-1,-25600,-1,-1

PIXELS 117,270 GET Chooser STYLE 65536 FONT "Chicago",12 PICTURE "6*RV Transcript profiles PIXELS 135,128 SAY Update STYLE 0 FONT "Geneva",12 SIZE 15,79 COLOR 0,0,0,-25600,-1,-1

PIXELS 171,128 SAY cloneno STYLE 0 FONT "Geneva",12 SIZE 15,79 COLOR 0,0,0,-25600,-1,-1

PIXELS 135,44 SAY "Last update: STYLE 65536 FONT "Geneva",12 COLOR 0,0,0,-1,-1,-1,-1,-1
@ PIXELS 133,44 SAY "Total clones: STYLE 65536 FONT 'Geneva',12 COLOR 0,0,-1,-1,-1,-1 @ PIXELS 45,296 SAY "V1.30" STYLE 65536 FONT 'Geneva',782 COLOR 0,0,-1,-1,-1,-1
 * EOF: Lifeseq menu.fmt
READ
DO CASE
CASE Chooser=1
DO "SmartGuy:FoxEASE+/Mac:fox files:Output programs:Master analysis 3.prg"
CASE Chooser=2
DO "SmartGuy:Fox3ASE+/Mac:fox files:Output programs:Subtraction 2.prg"
CASE Chooser=3
DO "SmartGuy:FoxBASE+/Mac:fox files:Output programs:Northern (single).prg"
CASE Chooser=4
USE *Libraries.dbf*
BROWSE
                                        12. 33.4
CASE Chooser=5
DO "SmartGuy:FoxEASE+/Mac:fox files:Output programs:See individual clone.prg"
DO "SmartGuy:FoxBASE+/Mac:fox files:Libraries:Output programs:Menu.prg"
CASE Chooser=7
CLEAR
SCREEN 1 OFF
RETURN
ENDCASE
LOOP
```

? 'Arranged by FUNCTION'

```
61,30 SAY "Database Subset Analysis" STYLE 65536 FONT "Geneva",274 COLOR 0,0,0,-1,-1,-1
30000 B
7<sub>KC</sub>(-5)
BANK BE IN A TO STORE
      "快点"、1967年
 ?? TIME()
? 'Clone numbers '
 ?? STR(INITIATE, 6, 0)
77 through '?? STR(TERMINATE, 6,0)
IF ENTIRE=1
7 All libraries
EMDIF
IF ENTIRE=2
MARK-1
DO WHILE .T.
  IP MARK>STOPIT
EXIT
ENDIF
USE SPLECTED
GO MARK
77 TRIM(libname)
    STORE MARK+1 TO MARK
    LOOP
? 'Designations: 'IF Ematch-0
    ENDDO
IF Ematch=0 .AND. Hmatch=0 .AND. Omatch=0
?? 'All'
               and the second of the second of
ENDIF
IF Ematch=1
?? 'Exact,'
                  ENDIF
IF Hmatch=1
?? "Human, '
ENDIF
IF Omatch=1
?? 'Other sp.'
ENDIF
IF CONDEN-1
? 'Condensed format analysis'
ENDIP
IF:ANAL=1
? 'Sorted by NUMBER'
ENDIF
IF ANAL=2
? 'Sorted by ENTRY'
ENDIF
IF ANAL=3
? 'Arranged by ABUNDANCE'
ENDIF
IF ANAL=4
? 'Sorted by INTEREST'
ENDIF
IF ANAL=5
? 'Arranged by LOCATION'
ENDIP
IF ANAL-5
? 'Arranged by DISTRIBUTION'
ENDIF
IF ANAL-7
```

```
PNDIF
? 'Total clones represented: '
?7 STR(STARTOT, 6, 0)
                                ? Total clones analyzed; ?? STR(ANALTOT, 6,0)
                                7 138 Type 10 10
                                                                                    3.00
                                                                           First Control of the 
                                                                      CBARTO DE
EXITADO EN LAS ACTUARDADAS MANAGO
PARTICIPATA DE CARROS ARA
MARCOLLAS ACTUARDADAS
MARCOLLAS ACTUARDADAS
SARTOS ACTUARDADAS
SARTOS ACTUARDADAS ACTUARDADAS
SARTOS ACTUARDADAS ACT
And the second of the second o
                                                                           8 100 CO.
                                                                                                                                                                                       And the second s
```

USE TEMP1 COUNT TO TOT COUNT TO TOT

?? 'Total of'

?? STR(TOT, 4:0)

?? 'Clones'

"list off fields number, L, D, F, Z, R, C, ENTRY, DESCRIPTOR, LENGTH, RFEND, INIT, I
list off fields number, L, D, F, Z, R, C, ENTRY, DESCRIPTOR

COORD Damage coor list off fields number, I CLOSE DATABASES ERASE TEMPI DBF USE TEMPDESIG

 $(W_{ij}) = (W_{ij}) \cup (W_{ij})$

8 1

The first that was a first time of the second

{ i

in the self to one your a mortalist control again of The face of the second management and second of the transfer of

the rever was a sine of the lar of the segmenters

att december of the a same of the attention got a sequence whose USE TEMP1 COUNT TO TOT COUNT TO TOT OF THE CONTROL AND ADDRESS OF THE PROPERTY OF THE CONTROL OF THE CON ?? STR(TOT, 4,0), transport that the body of the body section is a section of the ?? ' clones!

*list off fields number, L, D, F, Z, R, C, ENTRY, DESCRIPTOR, LENGTH, RFEND, INIT, I list off fields number, L, D, F, Z, R/C, EVTRY, DESCRIPTOR (1986) Sand Association CLOSE DATABASES

ERASE TEMPLOBER AND PROMISE TO A DESCRIPTION OF ANY AND ANY AN

USE TEMPOESIG Some the state of the control of the section of a first of the entropy of refer ere and at the fenerality sequences, there each erec in it is not the control of the body out to be say much

and the entering of the contract of the contract of the TO SEC TO A CONTRACTOR OF THE CONTRACTOR AND A SECURITION OF THE SECURITIES OF And the second section is the second section of the second

The second of the second of the second of the second of The desired the second of the second sections. and the second property of The second of the second of the second

Camarina Carlos and American State of the Carlos of the Ca , .

The first of the second second

A Commence of the Commence of

the Control of the second of t

. . .

en termination of a major of the second of t

```
*Northern (single), version 11-25-94
     Close, databases .

SEM MAIN APP TO A MILLION OF COMPLICATION OF SHORT OF COMPLETE SPECIAL 
     SET PRINT OFF
                                                                                  Commercial Section & Section Commence of the C
      SET EXACT OFF
     CLEAR
                                                                                TO Bobject have a think approximate accommon to the
     STORE '
    STORE 0 TO Numb
                                                                                                                                                                                                   ' TO Dobject
    STORE 0 TO Zog the apademistrate segment intergrate of printing part of storing property.
    * Program.: Northern (single).fmt
    * Date...: 8/ 8/94 Date...: FoxBASE+/Mac, revision 1.10
     * Notes. ... Format file Northern (single) and second confidence and
 SCREEN 1 TYPE 0 HEADING "Boreen 1" AT 40,2 SIZE 286,492 PIXELS FONT "Geneva",12 COLOR 0,0,0 PIXELS 15,81 TO 46,397 STYLE 28447 COLOR 0,0,-1,-25600,-1,-1 PIXELS 99,79 TO 192,422 STYLE 28447 COLOR 0,0,0,-25600,-1,-1 PIXELS 115,98 SAY "Entry #:" STYLE 65536 FONT "Geneva",12 COLOR 0,0,0,-1,-1,-1 PIXELS 115,173 GET Ecolect STYLE 0 FONT "Geneva",12 SIZE 15,142 COLOR 0,0,0,-1,-1,-1 PIXELS 145,89 SAY "Description" STYLE 65536 FONT "Geneva",12 COLOR 0,0,0,-1,-1,-1 PIXELS 145,173 GET Dobject STYLE 0 FONT "Geneva",12 SIZE 15,241 COLOR 0,0,0,-1,-1,-1 PIXELS 35,89 SAY "Single Northern search screen" STYLE 65536 FONT "Geneva",274 COLOR 0,0,-1,-1,-1 PIXELS 220,162 GET Bail STYLE 65536 FONT "Chicago",12 PICTURE "2*R Continue;Bail out SIZE PIXELS 175,98 SAY "Clone #:" STYLE 65536 FONT "Geneva":12 COLOR 0,0,0,-1,-1,-1
G PIXELS 175,98 SAY "Clone #:" STYLE 65536 FONT "Geneva";12 COLOR 0,0,0,-1,-1,-1
G PIXELS 175,173 GET Numb STYLE 0 FONT "Geneva",12 SIZE 15,70 COLOR 0,0,0,-1,-1,-1
G PIXELS 80,152 SAY "Enter any ONE of the following: STYLE 65536 FONT "Geneva",12 COLOR -1,
   * EOF: Northern (single).fmt
                                                                                                                                                                                                                                                             The water than the course of the water
  READ
                                                                                                          The first of the first of the state of the first of the state of the s
   IF Bail=2
   CLEAR
  screen 1 off
                                                                                                             Control to the control of the
                                                                                                                                                                                                                                                                                                                               Control of the second of the second
  USE "SmartGuy:FoxBASE+/Mac:Fox files:Lookup.dbf"
  SET TALK ON -
                                                                                                                                         A CARLO CONTRACTOR AND CONTRACTOR OF THE CONTRACTOR
  IF Eobject<>'
  STORE UPPER (Eobject) to Eobject
  SET SAFETY OFF
  SORT ON Entry TO "Lookup entry.dbf"
  SET SAFETY ON
  USE "Lookup entry.dbf"
 LOCATE FOR Look=Eobject
  IF .NOT.FOUND()
  CLEAR
 LOOP
 ENDIF
 BROWSE
  STORE Entry TO Searchval
 CLOSE DATABASES
ERASE "Lookin entry.dbf"
ENDIF
 IF Dobject >'
SET EXACT OFF
 SET SAFETY OFF
SORT ON descriptor TO "Lookup descriptor.dbf"
 SET SAFETY On
USB *Lookup descriptor.dbf*
LOCATE FOR UPPER(TRIM(descriptor))=UPPER(TRIM(Dobject))
IF .NOT.FOUND()
CLEAR
```

```
The state of the second second
                                     Therefore the control of the control of the state of the control o
                                                                                     TO THE WEST AND THE THE MEET WAS THEFT TO BE
                                   A Harry Commission will be highly governor give
                                                             to suplie the mean training him andones to determine
    LOOP
    ENDIF
                                                    CONTROL OF NOTICE OF THREE LEADS OF THE WALLS COL
    BROWSE
   STORE Entry TO Searchval
    CLOSE DATABASES
ERASE "Lookup descriptor.dbf"
    SET EXACT ON
    ENDIP
                                                                                        Liberarchic mater which remarks a property of
                                                                 Le de 151 imano
  USE *SmartGuy:FoxBASE+/Mac:Fox files:clones.dbf*
   GO Numb
                                                                                       ten Harri Committelia, of while amalem late late less
   BROWSE
    STORE Entry To Searchval
  ENDIP
                                                                                     to the law weeks, those enemy the higher were
  CLEAR
   ? 'Northern analysis for entry ' the control of the
   ?? Searchval
                                                                                   one to laber massers of all was consensually
  ? 'Enter Y to proceed'
                                                                                                                                                                  Park Style Bellet
  WAIT TO OK
  CLEAR
                                                                                     A TOP OF THE CONTROL THE INDIVIDUAL TO GREEN WAS
  IF UPPER (OK) <> 'Y'
  screen 1 off
                                                                                                                                                     CONTRACTOR STATE
  RETURN
  ENDIF
                                                                                                                       1 74 TV SEC.
                                                                                                     · 1.
                                                                                                                                                                                                      44.
                                                                                                                                                                                                                                                          THE SECTION OF STAR
  * COMPRESSION SUBROUTINE FOR Library, dbf
                                                                                                                                                                       4
                                                                                                                                                                                  Contraction of the Contraction
  ? 'Compressing the Libraries file now...'
 USE 'SmartGuy:FoxBASE+/Mac:Fox files:libraries.dbf'
  SET SAFETY OFF
  SORT ON library TO "Compressed libraries.dbf"
  * FOR entered>0
                                                                                       The second second section of the second second sections and the second
 SET GAFETY ON
 USE "Compressed libraries.dbf"
                                                                                                                                                                                      Charles to grant
 DELETE FOR entered=0
 PACK
 COUNT TO TOT
 MARK1 = 1
 8W2=0
 DO WHILE SW2=0 ROLL
        IF MARK1 >= TOT
        PACK
        SW2=1
       LOOP
                                                                                                                                                                                                                                                                                   . ...
       ENDIF
 GO MARKI
STORE library TO TESTA
 SKIP
 STORE Library TO TESTS
 IF TESTA = TESTB
DELETE
ENDIF
MARK1 = MARK1+1
 LOOP
ENDDO ROLL
* Northern analysis
CLEAR
? 'Doing the northern now...'
SET TALK ON
USE *SmartGuy:FoxBASE+/Mac:Fox files:clones.dbf*
SET SAFETY OFF
COPY TO "Hits.dbf" FOR entry=searchval
SET SAFETY ON
```

DO "Test print.prg"

RETURN

```
Which are plain all and marriages to 200 and about
        wants of one representing a programmer
CLOSE DATABASES (0) 1803 A LAC & DEPRESENTER LIVE POPULATION OF
SELECT 1
USE "Compressed libraries.dbf"
STORE RECOUNT(); TO Entries (10) Varies (20) 100 (20) (20) (20) (20) (20)
USE "Hits/Obf" payed nia or it, is sugar to magnific markets which
Mark=1
DO WHILE IT FOLLS ter gen from Joloh She bruche show the transcribed:
IF Marlo Entries (7) - Source among a common of Linear mach approves
EXIT
       time of lates within the take their of sleeps in me
ENDIF
GO MARK
STORE library TO Jigger of Parenting of the States
                                                  35.10
SELECT 2
COUNT TO Zog FOR library=Jigger Library=Land Country A we missing markets of educations of the
SELECT 1 (1) 20 Find whaten the stip production of the good in appearing the
Mark=Mark+1 :---
LOOP
ENDDO .
                   The house of the facilities
                                                   1828 B. G. GARDEN
SELECT 1
BROWSE FIELDS LIERARY, LIENAME, ENTERED, HITS AT 0,0
CLEAR
? 'Enter Y to print:'
                                     and the same
WAIT TO PRINSET
IF UPPER (PRINSET) = 'Y'
                             Part Para . . . . .
SET PRINT ON
CLEAR
                           Carthal Angline & Sanda Ar.
EJECT.
SCREEN 1 TYPE 0 HEADING "Screen 1" AT 40,2 SIZE 286,492 PIXELS FONT "Geneva",14 COLOR 0,0,0
? 'DATABASE ENTRIES MATCHING ENTRY '
?? Searchval
                                                                       Say Barret
? DATE()
                   All Carried
                               ...
SCREEN 1 TYPE 0 HEADING "Screen 1" AT 40,2 SIZE 286,492 PIXELS FONT "Geneva",7 COLOR 0,0,0,
LIST OFF FIELDS library, libname, entered, hits
SELECT 2
LIST OFF FIELDS NUMBER, LIBRARY, D.S.F.Z.R. ENTRY, DESCRIPTOR, RFSTART, START, RFEND
SET TALK OFF
SET PRINT OFF
ENDIP
CLOSE DATABASES
SET TALK OFF
CLEAR
```

The second one we will have a range the less was about the second with the sec

```
1996年,2016年,1986年,1986年,1986年,1986年,1986年,1986年,1986年,1986年,1986年,1986年,1986年,1986年,1986年,1986年,1986年,1986年,19
                                                                                                                                                           TABLE 6
                                                                                                                                                                                                                                                       CARLOS LANGER CARLO ARRIVANCE.
                                                To the company of the common content of the common that a content of the content 
                                             Charles and Charles and Charles and Charles
                                                                                                                                                                                        "我就不过有一个的结婚,这些女性,这一致好不少不多的好"而说我的这样说道。
                                            Normal in 1936 they are encountries are later for a page as
                                             of a tell balary a coffeeable the coloradio.
  library
                                            libname
                                            Inflamed adenoid and in the state of the sta
  ADENINB01
                                                                                                                                                                                                                                                                 Commercial Contraction
  ADRENOR01
                                            Adrenal gland (r)
                                                                                                                                                                                                              The grown also be given
  ADRENOTO1
                                            Adrenal gland (T) > 1. St. eng. (1.11) 13.5.5.5
                                                                                                                                                                                          u 1
 AMERIOTO: Access grams (1)

AMERIOTO: AMERICAN A
  CARDNOTO1 Cardiac muscle (T)
                                                                                                                  of a country of about 1 the efficient of a specimen
  CHAONOTO1 Chin. hamster overy
  CORNNOTO1 Corneal stroma FibraGT01 Fibroblest, AT 5
                                                                                                                                                                                                      1275 623
  FIBRAGT02
                                            Fibroblast, AT 30
 FIERANTO1
                                            Fibroblast, AT
                                            Fibroblast, uy 5
  FIERNGT01
                                                                                                                                                                            and other in some promiselying.
                                                                                                         44.0 By-1.0
 FIBRNGT02
                                            Fibroblast, uv 30
                                            Fibroblast . .
  FIRRNOTO:
                                            Fibroblast, normal
                                                                                                                                                                                                     100 miles 188 188
                                                                                                                                                                                                                                                                     The Article Specification of
 FIZANOTO2
HMC1NOT01 Mast cell line HMC-1
HUVELPB01 HUVEC IFN,TNF,LPS
                                                                                                                                                                                the country like one cases of an execution
                                                                                                                         CONTRACT
  HUVENOBO1 HUVEC control
 HUVESTB01 HUVEC shear stress
  HYPONOB01 Hypothelamus
 KIDNNOT01 Kidney (T)
                                            Liver (T)
LIVRNOT01 Liver (T)
LUNGNOT01 Lung (T)
                                                                                                                              1.
                                                                                                                                                                                                          Section 1
 MUSCNOT01 Skaletal muscle (T)
 OVIDNOBO1 Oviduct
 PANCNOTO: Pancreas, normal
 PITUNOR01
                                            Pituitary (r)
 PITUNOT01
                                            Pitultary (T)
 PLACNOBO1 Piscenta
  SINTNOTO2
                                            Small intestine (T)
                                                                                                                   ~,.
 SPLNFET01
                                            Spleen+liver, fetal
SPLNNOTO2 Spleen (T)
STOMNOTO1 Slomach
 6YNORAE01 Rheum synovium
 TELYNOTO1
                                            T + B lymphoblast
                                                                                                                                                                                                                                                                                                                                                                              Same State
TESTNOTO:
                                            Testis (T)
THP1NOB01
                                          THP-1 control
THP1PEB01
                                           THP phorbol
THP1PLR01
                                           THP-1 phorbol LPS
UB37NOT01 US37, monocytic leuk
```

numbe	rlibrary	d	2	f	2 7	entry	descriptor	rfetari	81871	rfend
2304	U937NOT01	E	н	C	CT	HUMEF#B	Elongation lactor 1-beta	Ď.	0	773
3240	HMC1NOT01	E	Н	C	C T	HUMEF1B	Elongation factor 1-beta	ŏ	370	773
3259	HMC1NOT01	E	н	C	C T	HUMEF18	Elongation factor 1-beta	ō	371	773
4693	HMC1NOT01	E	Н	C	C T	HUMEF1B	Elongation factor 1-beta	ō	470	773
8989	HMC1NOT01	Ε	н	C	C T	HUMEF1B	Elongation factor 1-beta	ō	327	773
9139	HMC1NOT01	Ε	H	C	CT	HUMEF18	Elongation factor 1-beta	٥	375	773

WHAT IS CLAIMED IS:

1. A method of analyzing a specimen containing gene transcripts, said method comprising the steps of:

- (a) producing a library of biological sequences;
- (b) generating a set of transcript sequences, where each of the transcript sequences in said set is indicative of a different one of the biological sequences of the library;
- programmed computer in which a database of reference transcript sequences indicative of reference biological sequences is stored, to generate an identified sequence value for each of the transcript sequences, where each said identified sequence value is indicative of a sequence annotation and a degree of match between one of the transcript sequences and at least one of the reference transcript sequences; and
- (d) processing each said identified sequence value to generate final data values indicative of a number of times20 each identified sequence value is present in the library.
 - 2. The method of claim 1, wherein step (a) includes the steps of:

obtaining a mixture of mRNA; making cDNA copies of the mRNA;

- isolating a representative population of clones transfected with the cDNA and producing therefrom the library of biological sequences.
 - 3. The method of claim 1, wherein the biological sequences are cDNA sequences.
- 4. The method of claim 1, wherein the biological sequences are RNA sequences.
 - 5. The method of claim 1, wherein the biological sequences are protein sequences.

6. The method of claim 1, wherein a first value of said degree of match is indicative of an exact match, and a second value of said degree of match is indicative of a non-exact match.

- 7. A method of comparing two specimens containing gene transcripts, said method comprising:
 - (a) analyzing a first specimen according to the method of claim 1;
- (b) producing a second library of biological10 sequences;
 - (c) generating a second set of transcript sequences, where each of the transcript sequences in said second set is indicative of a different one of the biological sequences of the second library;
- (d) processing the second set of transcript sequences in said programmed computer to generate a second set of identified sequence values known as further identified sequence values, where each of the further identified sequence values is indicative of a sequence annotation and a degree of match between one of the biological sequences of the second library and at least one of the reference sequences;
- (e) processing each said further identified sequence value to generate further final data values indicative of a
 25 number of times each further identified sequence value is present in the second library; and
- (f) processing the final data values from the first specimen and the further identified sequence values from the second specimen to generate ratios of transcript
 sequences, each of said ratio values indicative of differences in numbers of gene transcripts between the two specimens.
- 8. A method of quantifying relative abundance of mRNA in a biological specimen, said method comprising the steps of:
 - (a) isolating a population of mRNA transcripts from the biological specimen;

(b) identifying genes from which the mRNA was transcribed by a sequence-specific method;

- (c) determining numbers of mRNA transcripts corresponding to each of the genes; and
- (d) using the mRNA transcript numbers to determine the relative abundance of mRNA transcripts within the population of mRNA transcripts.
 - 9. A diagnostic method which comprises producing a gene transcript image, said method comprising the steps of:
- (a) isolating a population of mRNA transcripts from a biological specimen;
 - (b) identifying genes from which the mRNA was transcribed by a sequence-specific method;
- (d) using the mRNA transcript numbers to determine the relative abundance of mRNA transcripts within the population of mRNA transcripts, where data determining the relative abundance values of mRNA transcripts is the gene transcript image of the biological specimen.
 - 10. The method of claim 9, further comprising:
 - (e) providing a set of standard normal and diseased gene transcript images; and
- (f) comparing the gene transcript image of the biological specimen with the gene transcript images of step (e) to identify at least one of the standard gene transcript images which most closely approximate the gene transcript image of the biological specimen.
- 11. The method of claim 9, wherein the biological specimen is biopsy tissue, sputum, blood or urine.
 - 12. A method of producing a gene transcript image, said method comprising the steps of
 - (a) obtaining a mixture of mRNA;
 - (b) making cDNA copies of the mRNA;

(c) vinserting the cDNA into a suitable vector and using said vector to transfect suitable host strain cells which are plated out and permitted to grow into clones, each clone representing a unique mRNA;

- 5 (d) isolating a representative population of recombinant clones;
 - (e) identifying amplified cDNAs from each clone in the population by a sequence-specific method which identifies gene from which the unique mRNA was transcribed;
- (f) determining a number of times each gene is represented within the population of clones as an indication of relative abundance; and
- (g) listing the genes and their relative abundance in order of abundance, thereby producing the gene transcript
 15 'image.
 - 13. The method of claim 12, also including the step of diagnosing disease by:

repeating steps (a) through (g) on biological specimens from random sample of normal and diseased humans, encompassing a variety of diseases, to produce reference sets of normal and diseased gene transcript images;

obtaining a test specimen from a human, and producing a test gene transcript image by performing steps (a) through (g) on said test specimen;

comparing the test gene transcript image with the reference sets of gene transcript images; and

identifying at least one of the reference gene transcript images which most closely approximates the test gene transcript image.

30 14. A computer system for analyzing a library of biological sequences, said system including:

means for receiving a set of transcript sequences, where each of the transcript sequences is indicative of a different one of the biological sequences of the library;

35 and

means for processing the transcript sequences in the computer system in which a database of reference transcript

sequences indicative of reference biological sequences is stored, wherein the computer is programmed with software for generating an identified sequence value for each of the transcript sequences, where each said identified sequence value is indicative of a sequence annotation and a degree of match between a different one of the biological sequences of the library and at least one of the reference transcript sequences, and for processing each said identified sequence value to generate final data values indicative of a number of times each identified sequence value is present in the library.

15. The system of claim 14, also including:
library generation means for producing the library of
biological sequences and generating said set of transcript
15 sequences from said library.

16. The system of claim 15, wherein the library generation means includes:

静されられ さんだって Links Control Control Part point in the Control

20

means for obtaining a mixture of mRNA;

means for making cDNA copies of the mRNA;

means for inserting the cDNA copies into cells and
permitting the cells to grow into clones;

means for isolating a representative population of the clones and producing therefrom the library of biological sequences.

SYBASE database Structure

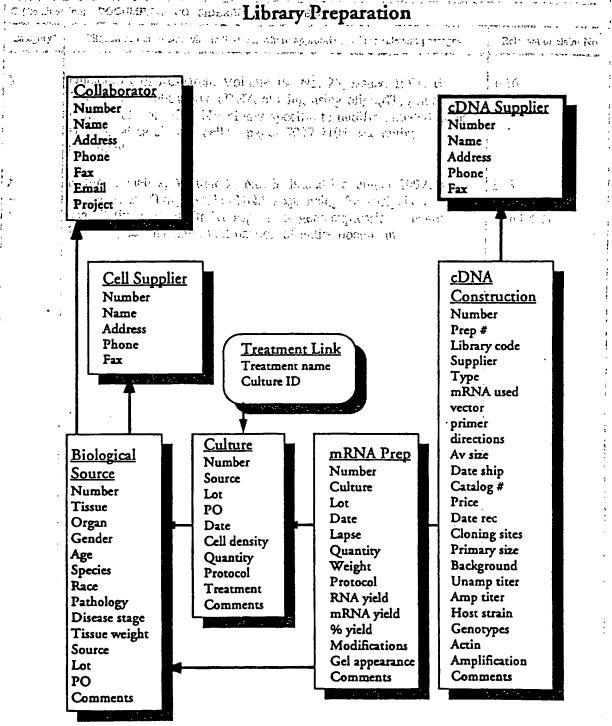


Figure 1

I Dentified Sequences Tempnum Tempred Templesig Templib Libsort Subsost Tempsubsost Temptorsort Final Data

Figure 2

WO 95/20681

PCT/US95/01160

Figure 3

The second of th

Incyte Bioinformatics Process

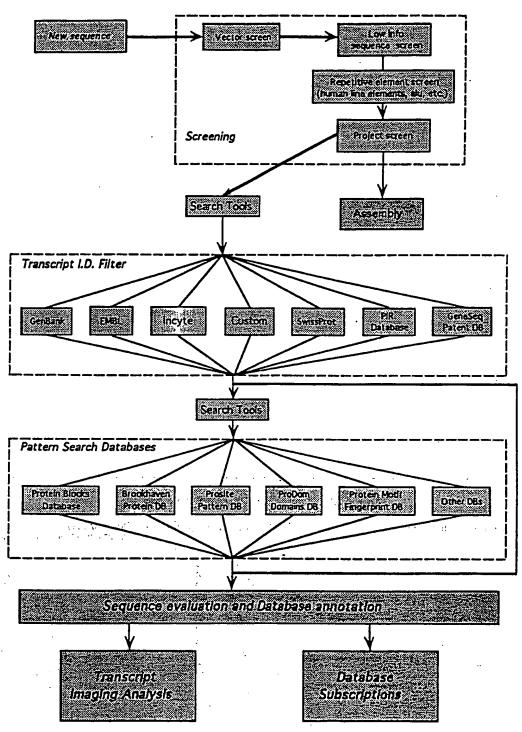


Figure 4

INTERNATIONAL SEARCH REPORT

International application No.
PCT/US95/01160

IPC(6) US CL	SSIFICATION OF SUBJECT MATTER : : : : : : : : : : : : : : : : : : :		· .					
								
		by classification symbols)						
U.S. :	ocumentation searched (classification system followed 435/6; 364/413.02	•						
٠.	The state of the s							
Documentat	tion searched other than minimum documentation to the	extent that such documents are included	in the fields searched					
	Professional States							
Electronic d	ata base consulted during the international search (nar	ne of data base and, where practicable	, search terms used)					
	LINE, APS, transcript, transcripts, cdan#; mrna#							
	more than well and leave the	was the tong of head	· Gor Dynten					
C. DOC		3. 1 3. 43. 15 (1.15) (1.15)						
Category*	Citation of document, with indication, where app	propriate, of the relevant passages	Relevant to claim No.					
X Y	IntelliGenetics Suite, Release 5.4, A issued January 1993 by IntelliGe Camino Real, Mountain View, C States of America, pages (1-6)-(1-entire document.	dvanced Training Manual, netics, Inc. 700 East El alifornia 94040, United 19) and (2-9)-(2-14), see	The state of the s					
?*	the second of the second of the second	pages 1651-1656, see	is at the particular con- months with the same					
ľ		Company of the state of the sta						
İ	and the second of the second of the second of	i	i					
			100000					
.	A Committee of the Comm							
								
X Furth	er documents are listed in the continuation of Box C.	See patent family annex.						
A* doc	cial categories of cited documents: ument defining the general state of the art which is not considered to of particular relevance	T leter document published after the inte date and not in conflict with the applica principle or theory underlying the inve	tion but cited to understand the					
	and the second s							
cite	ument which may throw doubts on priority claim(s) or which is d to establish the publication date of another citation or other	when the document is taken alone 'Y' document of particular relevance: the	. ala:al :al					
	cini reason (as specifical) umont referring to an oral disclosure, use, exhibition or other	'Y" document of particular relevance; the considered to involve an inventive combined with one or more other such	step when the document is					
mot	means being obvious to a person skilled in the art							
the	priority date claimed	'&' document member of the same patent						
27 APRIL	· ·	Date of mailing of the international sea 0 4 MAY 1995	rch report					
Commission Box PCT	per of Patents and Trademarks	JAMES MARTINELL Telephone No. (703) 308-0196	w trong					

Form PCT/ISA/210 (second sheet)(July 1992)*

INTERNATIONAL SEARCH REPORT

International application No. And the state of t

Category*	Citation of document, with indication, where appropri	ate, of the relevant passages	Relevant to claim No.
Y	Nucleic Acids Research, Volume 19, No. 2 Hara et al, "Subtractive cDNA cloning usin PCR: isolation of cDNA clones specific to embryonal carcinoma cells", pages 7097-71	5, issued 1991, E. g oligo(dT) ₃₀ -latex and undifferentiated human 04, see entire	1-16
	document.		
X Y	Nature Genetics, Volume 2, No. 3, issued of Okubo et al, "Large scale cDNA sequencing quantitative and qualitative aspects of gene (173-179, see narrative text portion of entire)	g for analysis of expression", pages	2 and 4-16
	on the transmission moderation of		
1 -	one e milio i dio milianga midal		
	the transfer of the transfer that the second		Company Car
	and the source of the second and second	A company of the contract of t	bojo prest
	te was a State of	mar Sirary Burn	0.15
	the second Attracts of the leading of the second	· ·	
,	at the with a compatible of	Comment of the second	. seine
		www.logosta.com	e e e e e e e e e e e e e e e e e e e
	·	The Contract	
İ	Description of the Control of the Artist	•	`
	et produktiva i saman katalog saman katalog saman katalog saman katalog saman katalog saman katalog saman kata Katalog saman katalog sama		for a figure
	 And the first of the control of the co		
I	to this community will be the constant of the second of the constant of the co	. •	er en en andere
	error ett er ett i skriver i skriver i skriver i skriver i skriver i skriver i skriver i skriver i skriver i s De skriver i skriver i skriver i skriver i skriver i skriver i skriver i skriver i skriver i skriver i skriver	•	
I	·		1 (1)
		· ·	
			1
	en de la companya de la companya de la companya de la companya de la companya de la companya de la companya de La companya de la co		ŧ
	* 44 *	· · · · · · · · · · · · · · · · · · ·	
	•	•	,
	$\ddot{\cdot}$		
		•	